



⁶⁸Ga-DOTATATE QC Analysis by Radiometric Detection

ACE[®]
Ultra-inert
UHPLC & HPLC Columns

Application #AN2720

Conditions

Column: ACE 3 C18
Dimensions: 150 x 3.0 mm
Part Number: ACE-111-1503
Mobile Phase: 0.1% TFA in H₂O/MeCN (76:24 v/v)
Flow Rate: 0.6 mL/min
Injection: 20 µL
Detection: Radiometric

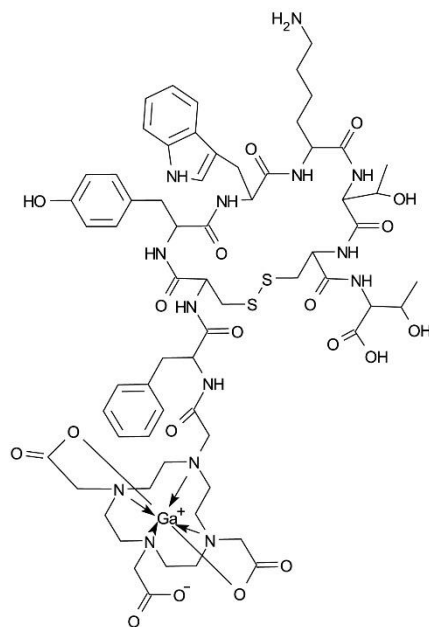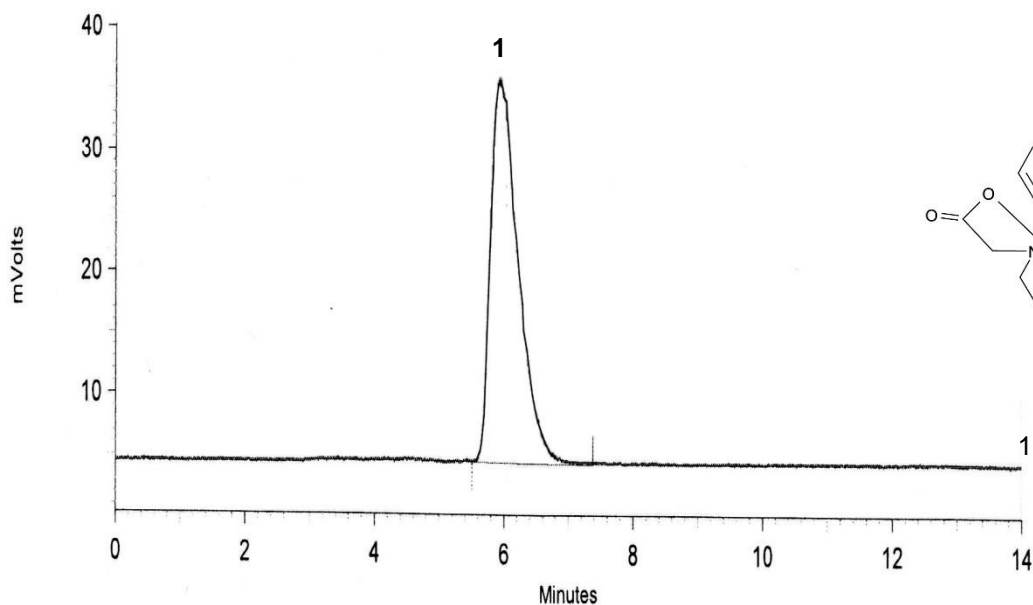

1. ⁶⁸Ga-DOTATATE

Reproduced with permission of Centre of Radiopharmaceutical Science, University Hospital Zurich, Switzerland

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com

ACE[®]
UHPLC & HPLC
Columns