


Steroids

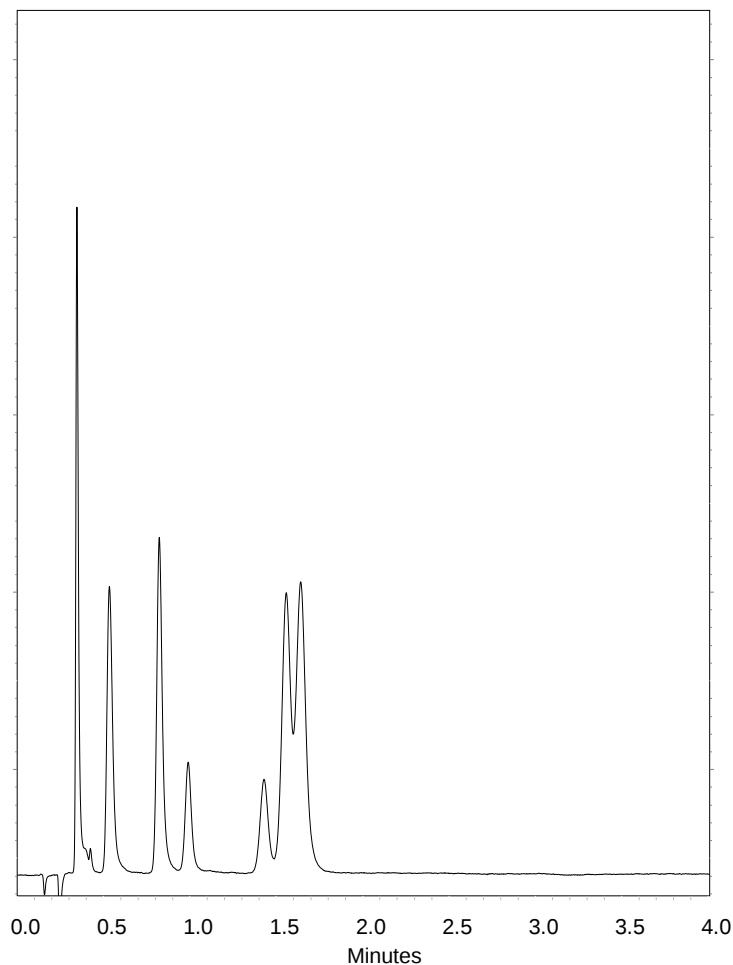

Column : **3 μ Fortis Cyano
50x2.1mm**

p/n : **FCN-020303**

Mobile Phase: **90:5:5 Heptane:THF:ACN**

Flow : **0.8ml/min**

Temp : **25°C**

λ : **254nm**

- 1. Progesterone**
- 2. 17 α Hydroxyprogesterone**
- 3. 11 α Hydroxyprogesterone**
- 4. Hydroxy-21-acetate**
- 5. Cortisone**
- 6. Prednisone**
- 7. Prednisolone**