


Determination of Alanyl Glutamine

ACE[®]
Ultra-inert
UHPLC & HPLC Columns


Conditions

Columns: ACE Excel 5 NH₂
Dimensions: 250 x 4.6 mm
Part Numbers: EXL-1214-2546U
Mobile Phase: 50 mM KH₂PO₄ in H₂O (pH 4)/MeCN (35:65 v/v)
Flow Rate: 0.7 mL/min
Injection: 20 µL
Temperature: 30 °C
Detection: UV, 215 nm
Method: ChP.2015, Vol. 2


Application #AN6570

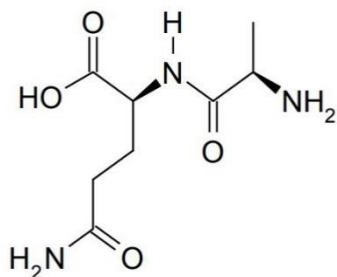

1. D-Alanyl-L-Glutamine

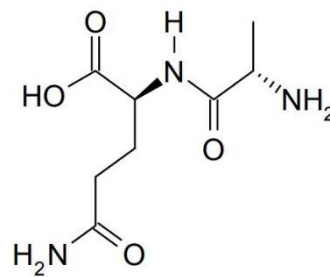


2. L-Alanyl-L-Glutamine

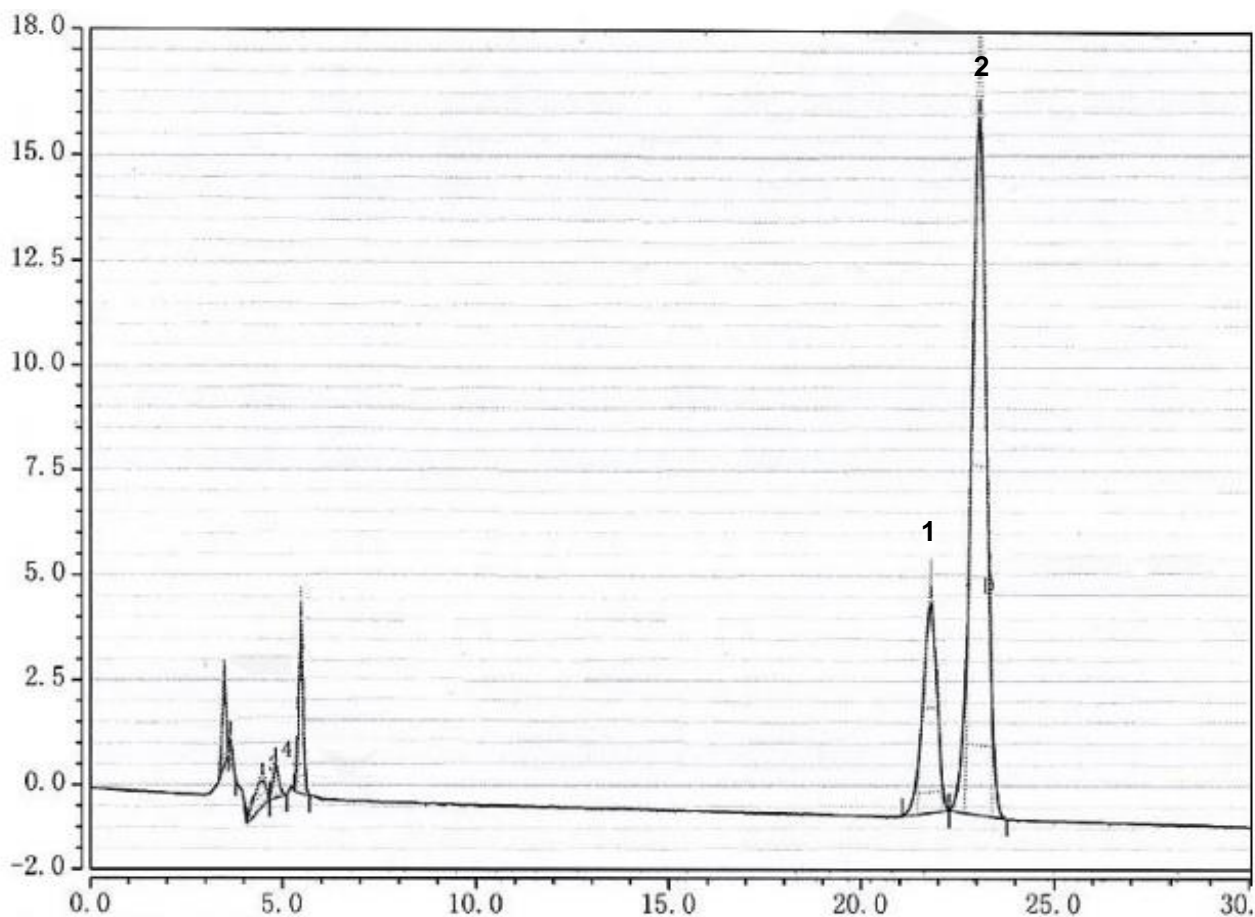

© 2018 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com


ACE[®]
UHPLC & HPLC
Columns