


## Application note #7630

# Separation of Alpha-, Beta- and Gamma-Cyclodextrin

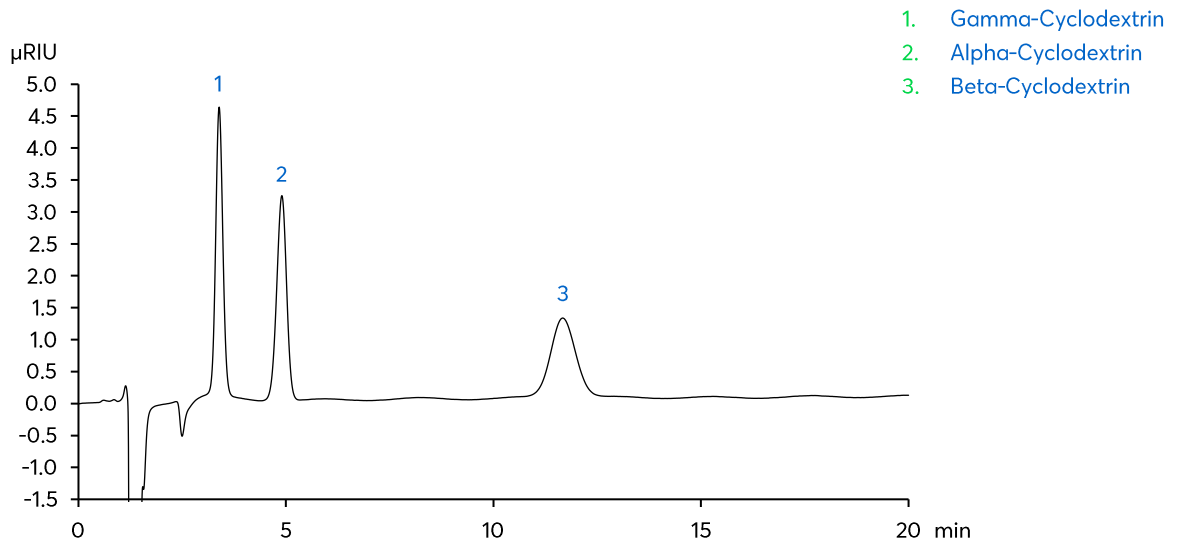

# Method Details

## CONDITIONS

Columns:	Avantor® ACE® Excel® 5 C18-Amide
Particle Size:	5 µm
Dimensions:	150 x 4.6 mm
Mobile Phase:	Water: Methanol (93:7)
Flow Rate:	1.5 mL/min
Injection Volume:	100 µL
Temperature:	30 °C
Detection:	RI

## ORDERING TABLE

Product	Details	Size	Part number
Avantor® ACE® Excel® 5 C18-Amide	HPLC Column	150 x 4.6 mm	EXL-1212-1546U

Reproduced with permission of Mumbai Avantor Application lab.