



Amino Acids

- Shimadzu Nexera
- A: 5mm ammonium acetate
- B: ACN

Shimpack VP-ODS
150x2.0mm 2.2um

1.0ml/min	
Time	%B
0	4
0.5	4
7.0	33
8.0	60
10	4

- 40oC
- 254nm

1.7um Fortis C18
50x2.1mm
30x2.1mm


0.9ml/min	
Time	%B
0	4
0.2	4
2.5	33
3.0	60
5	4

- 1) ASP
- 2) GLU
- 3) SER
- 4) GLY
- 5) HIS
- 6) ARG
- 7) THR
- 8) ALA
- 9) PRO
- 10) TYR
- 11) VAL
- 12) MET
- 13) CYS-CYS
- 14) ILE
- 15) LEU
- 16) PHE
- 17) LYS




Amino Acids

Shimpack VP-ODS
150x2.0mm 2.2µm

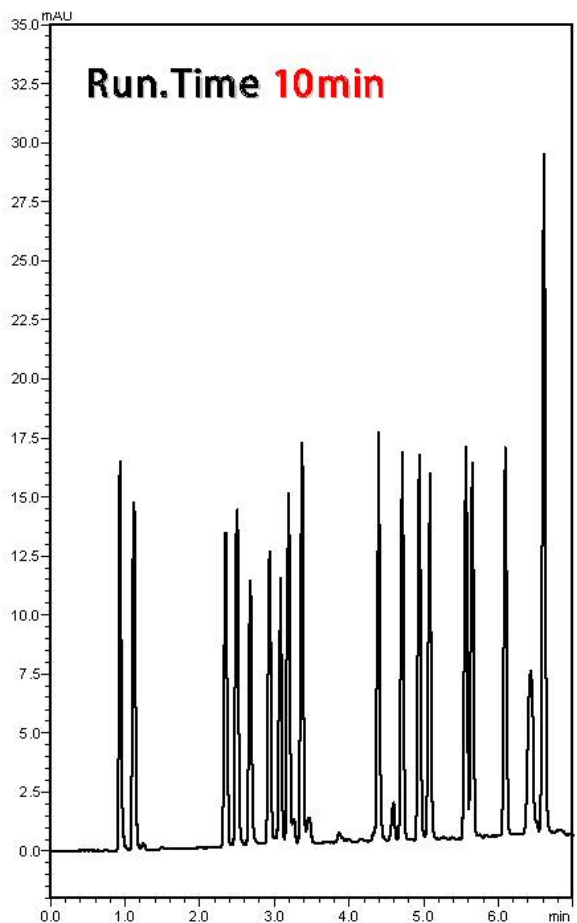

1.7µm Fortis C18
50x2.1mm

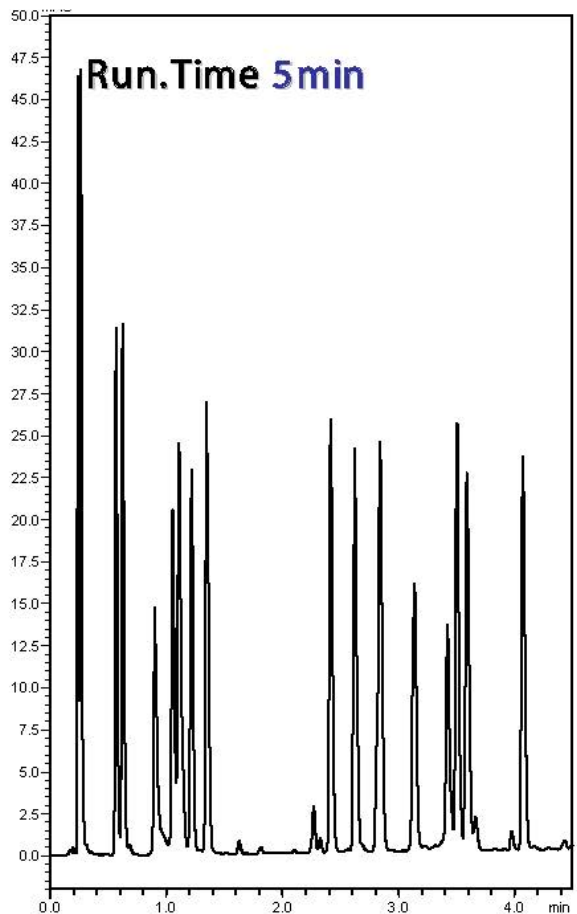

1.7µm Fortis C18
30x2.1mm

