

Amoxicillin and Clavulanate Potassium Tablets Related Substances

ACE[®]
Ultra-inert
UHPLC & HPLC Columns



Application #AN6270

Conditions

Column: ACE Excel 5 C18
Dimensions: 150 x 4.6 mm
Part Number: EXL-121-1546U
Mobile Phase: A: 10 mM KH₂PO₄ in H₂O (Adjusted to pH 6.0 with 2M NaOH)
B: 10 mM KH₂PO₄ in H₂O (Adjusted to pH 6.0 with 2M NaOH)/ MeCN (20:80)

Time (mins)	%B
0	2
20	30
22	2
32	2

Flow Rate: 1.0 mL/min
Injection: 20 µL
Temperature: 40 °C
Detection: UV, 230 nm
Monograph Method: ChP. 2015, Vol. 2

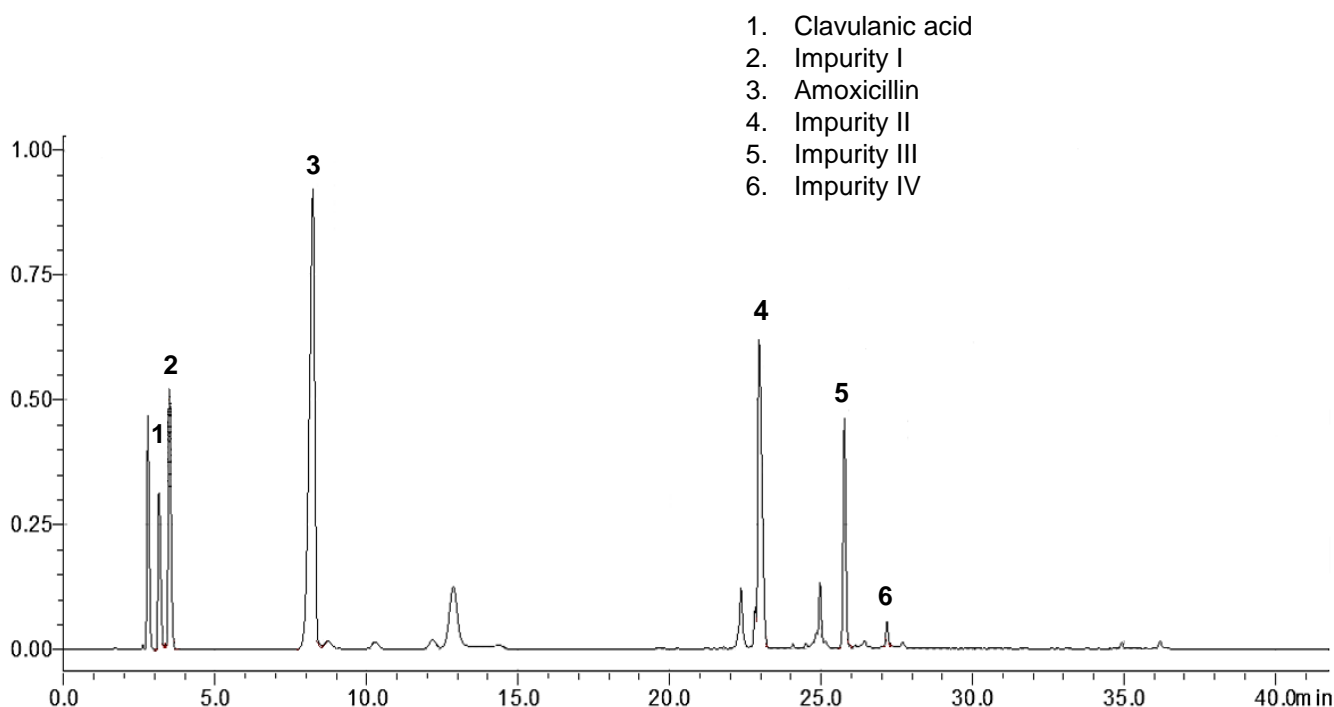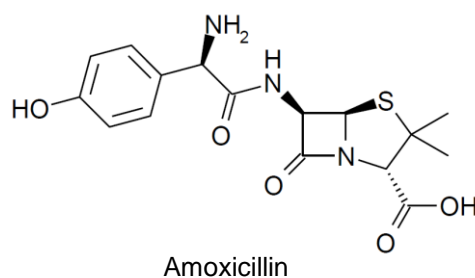

© 2018 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com



ACE[®]
UHPLC & HPLC
Columns