



Conditions

Column: ACE Excel 2 SuperC18
Dimensions: 50 x 2.1 mm
Part Number: EXL-1011-0502U
Mobile Phase: A: 0.1% formic acid in H₂O
B: 0.1% formic acid in MeCN


Time (mins)	%B
0.0	5
1.0	100
1.5	100
1.6	5

Post time 1.5 minutes


Flow Rate: 1.2 mL/min
Injection: 0.5 µL
Temperature: 50 °C
Detection: UV, 214 nm

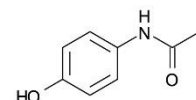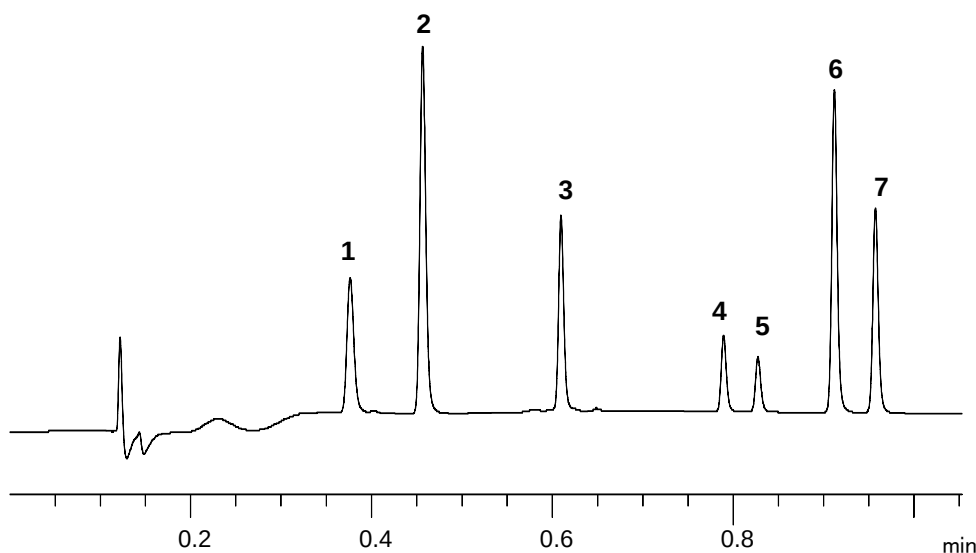


1. Paracetamol

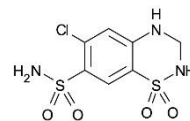


2. Hydrochlorothiazide

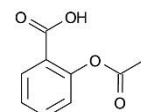


3. Aspirin

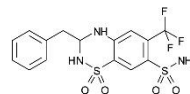


4. Bendroflumethiazide

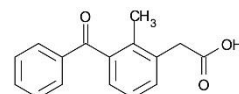

5. Ketoprofen

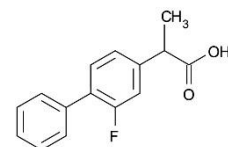

6. Flurbiprofen

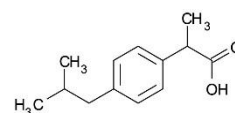

7. Ibuprofen