



Application #AN2150

Conditions


Column: ACE 5 C18-300
Dimensions: 250 x 4.6 mm
Part Number: ACE-221-2546
Mobile Phase: A: 0.1 % TFA in H₂O
B: 0.1% TFA in H₂O/MeCN (80:20 v/v)

Time (mins)	%B
0	25
15	40


Flow Rate: 1 mL/min
Temperature: Ambient
Detection: UV, 215 nm

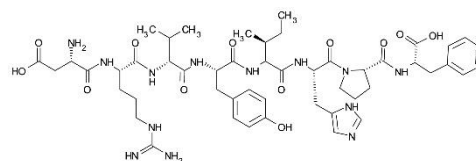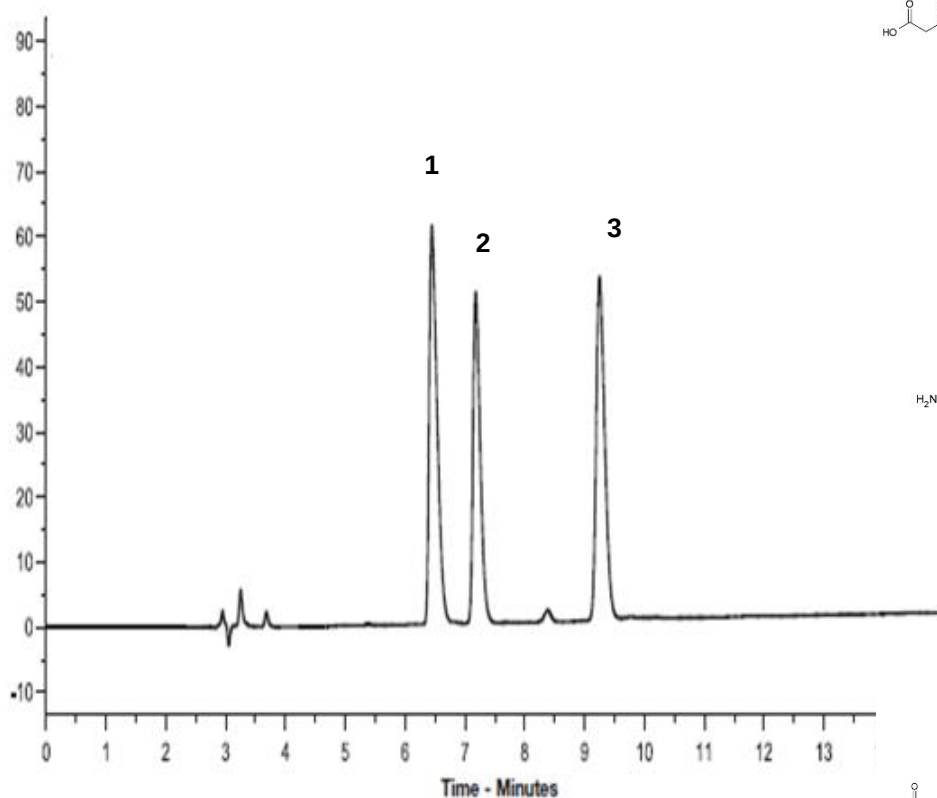

1. Angiotensin II
(*m/z* 1046.2)

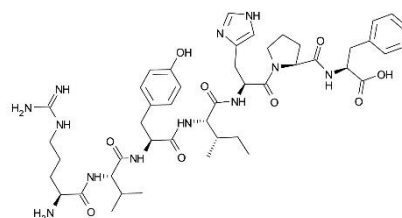

2. Angiotensin III
(*m/z* 931.1)

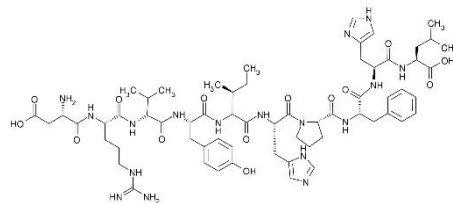

3. Angiotensin I
(*m/z* 1296.5)