


Conditions

Column: ACE 5 C18
Dimensions: 250 x 4.6 mm
Part Number: ACE-121-2546
Mobile Phase: 50 mM KH₂PO₄ pH 3.0/MeOH (50:50 v/v)
Flow Rate: 1 mL/min
Temperature: 22 °C
Detection: UV, 220 nm

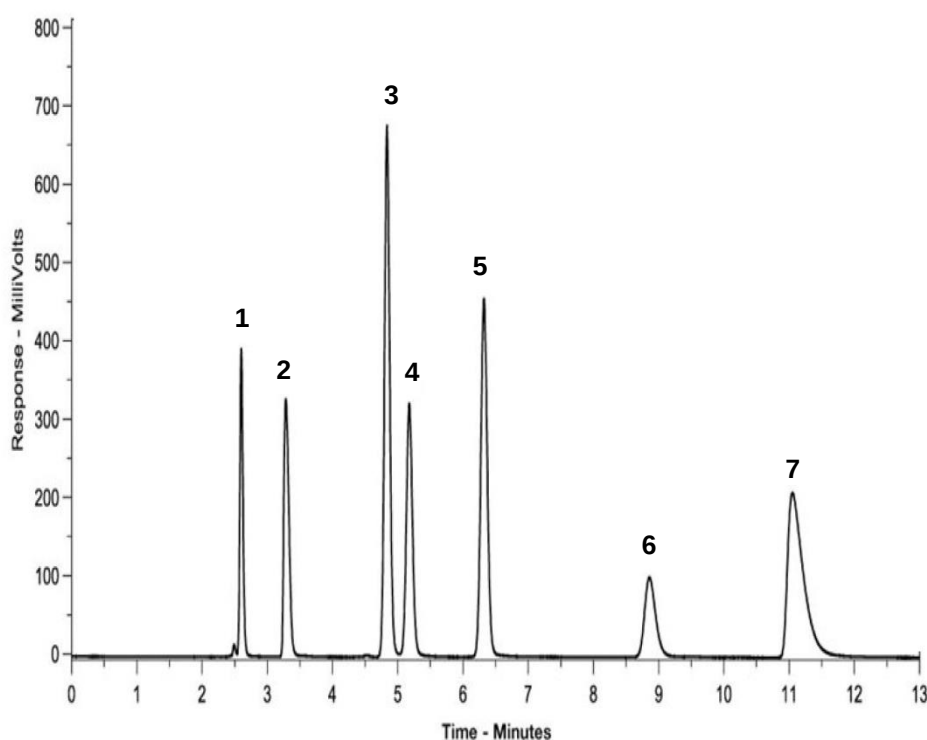

1. Maleic acid
2. Norephedrine
3. Salicylamide
4. Guaifenesin
5. Guaiacol
6. Chlorpheniramine maleate
7. Dextromethorphan