


BP Monograph - Lactulose Solution and Related Substances

ACE[®]
Ultra-inert
UHPLC & HPLC Columns

Application #AN4860

Conditions

Column: ACE Excel 3 NH₂
Dimensions: 250 x 4.6 mm
Part Number: EXL-1114-2546U
Mobile Phase: 0.253g sodium dihydrogen phosphate in 200 mL H₂O and diluted to 1000 mL with MeCN
Flow Rate: 1.1 mL/min
Injection: 20 µL
Temperature: 30 °C
Detection: Refractive Index

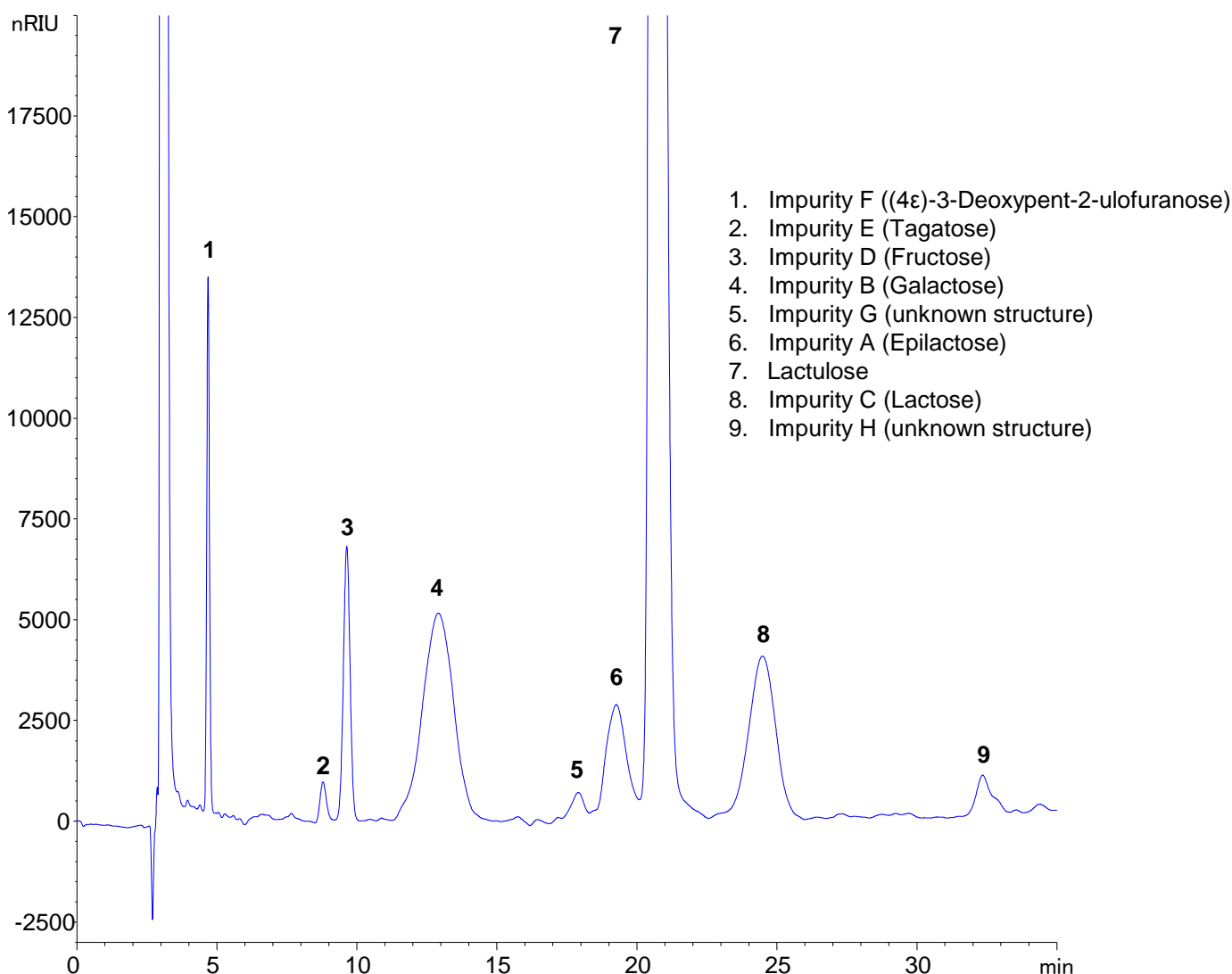

Reproduced with permission of Sterling Pharmaceuticals Ltd, Birmingham, UK

© 2017 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com

