


Conditions

Column: ACE 3 C18
Dimensions: 100 x 2.1 mm
Part Number: EXL-111-1002
Mobile Phase: A: 20 mM KH₂PO₄, pH 2.7
B: 20 mM KH₂PO₄, pH 2.7 in MeCN/H₂O (65:35 v/v)

Time (mins)	%B
0	5
1	5
12	95
13	95
14	5

Post time 3 minutes

Flow Rate: 0.6 mL/min
Injection: 2 µL
Temperature: 36 °C
Detection: UV, 214 nm

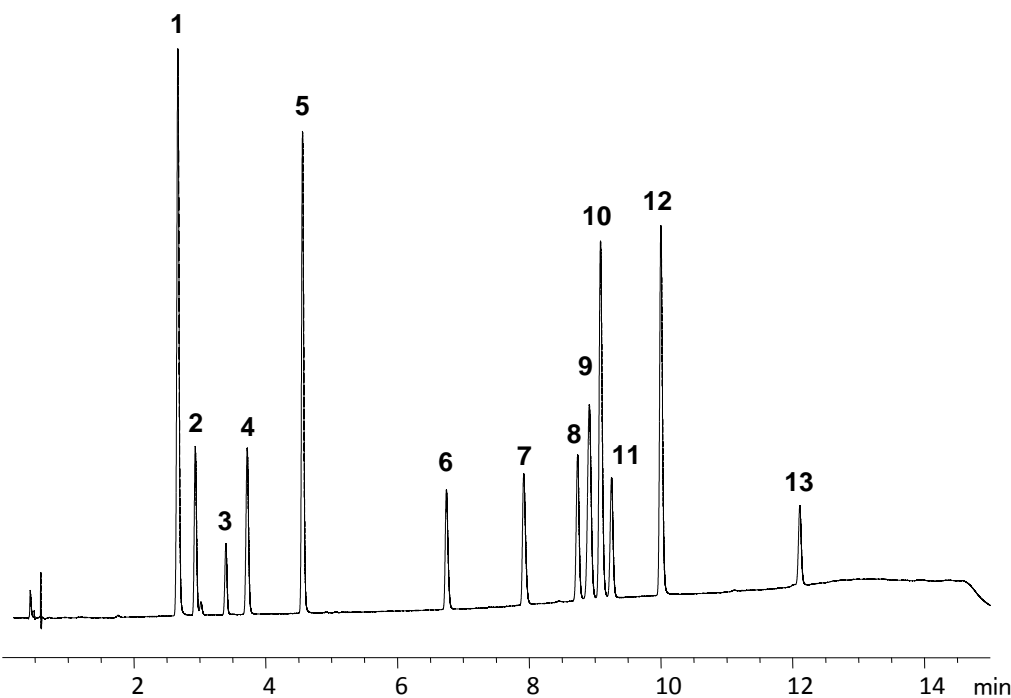

1. Hydrochlorothiazide
2. Sotalol
3. Amiloride
4. Atenolol
5. Pindolol
6. Metoprolol
7. Oxprenolol
8. Furesomide
9. Indapamide
10. Propranolol
11. Bendroflumethiazide
12. Carvedilol
13. Spironolactone