


Separation of β -thymidine (dT) and α -thymidine

ACE[®]
Ultra-inert
UHPLC & HPLC Columns


Application #AN6000

Conditions


Column: ACE Excel 5 C18-AR
Dimensions: 150 x 4.6 mm
Part Number: EXL-129-1546U
Mobile Phase: 10 mM ammonium acetate (aq.)
Flow Rate: 0.5 mL/min
Injection: 5 μ L
Temperature: 60 $^{\circ}$ C
Detection: UV, 260 nm

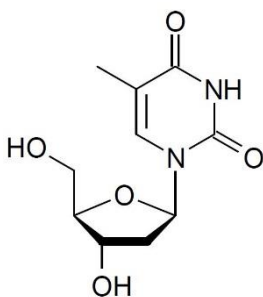

1. β -thymidine

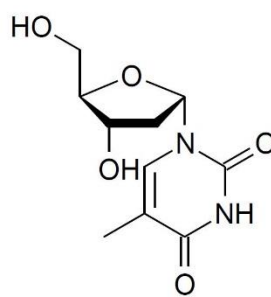

2. α -thymidine

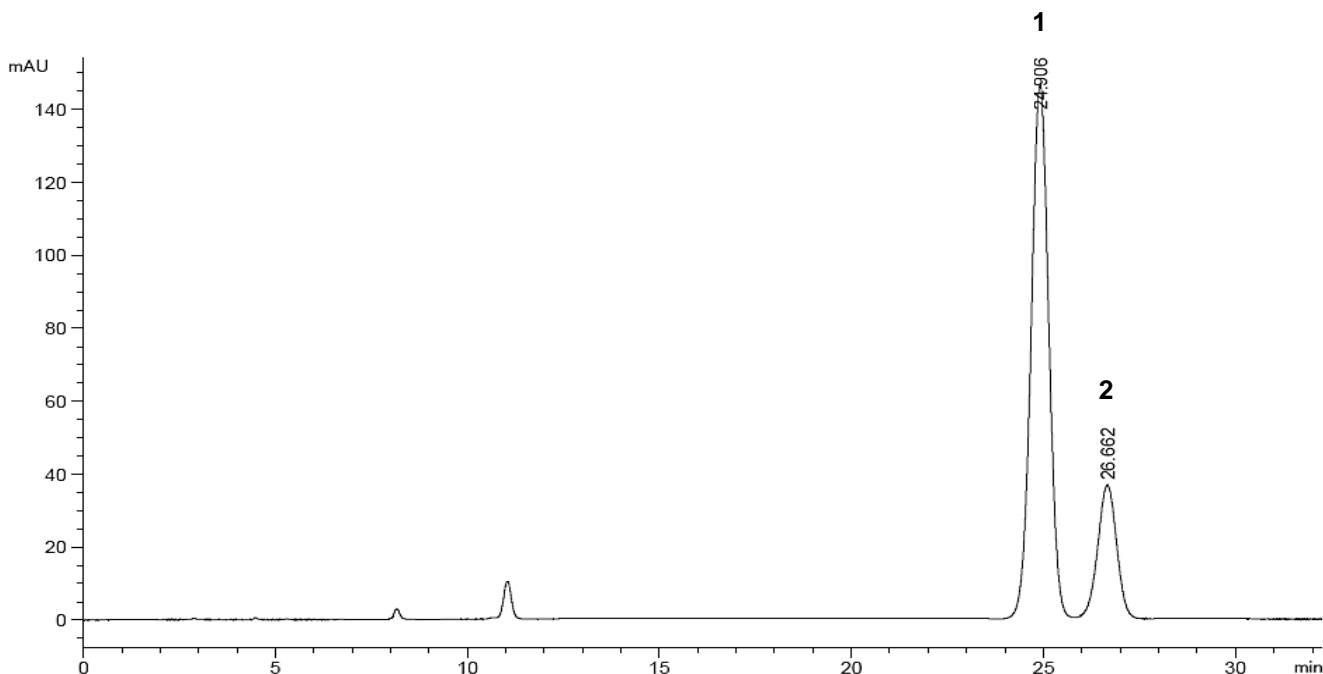

© 2018 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com


ACE[®]
UHPLC & HPLC
Columns