


Conditions

Column:	ACE 5 AQ
Dimensions:	150 x 4.6 mm
Part Number:	ACE-126-1546
Mobile Phase:	50 mM KH ₂ PO ₄ pH 3.0 in H ₂ O
Flow Rate:	1 mL/min
Injection:	2 µL
Temperature:	Ambient
Detection:	UV, 210 nm

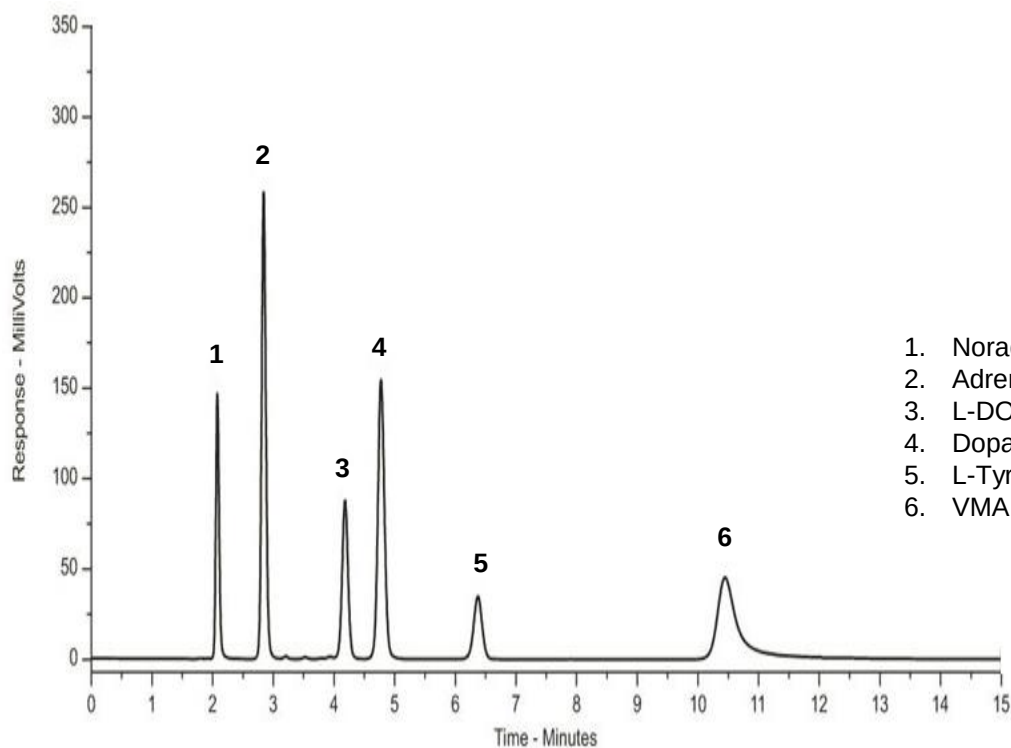

1. Noradrenaline
2. Adrenaline
3. L-DOPA
4. Dopamine
5. L-Tyrosine
6. VMA (Vanillylmandelic acid)