



Catecholamines and Metanephries Separation (Isocratic)

ACE[®]
Ultra-inert
UHPLC & HPLC Columns


Application #AN1490

Conditions


Column: ACE 5 C18-PFP
Dimensions: 150 x 4.6 mm
Part Number: ACE-1210-1546
Mobile Phase: 20 mM ammonium formate pH 3.0 in H₂O
Flow Rate: 1 mL/min
Injection: 5 µL
Temperature: 25 °C
Detection: UV, 260 nm

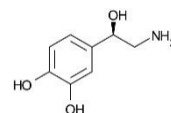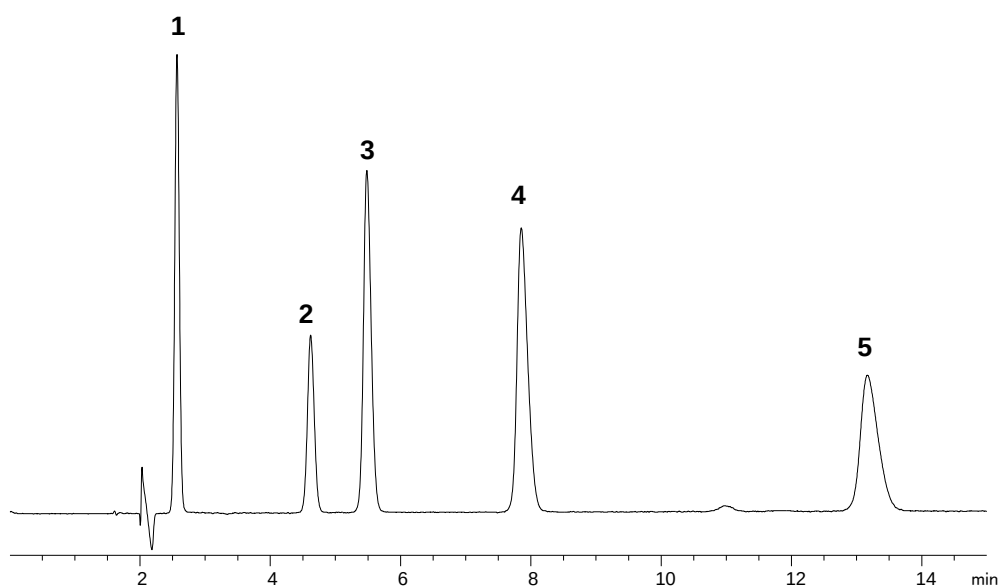


1. Norepinephrine

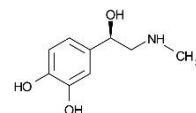


2. Epinephrine

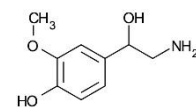

3. Normetanephrine

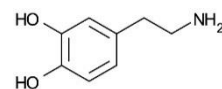

4. Dopamine

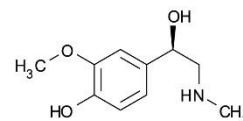

5. Metanephrine