


Clausena Excavate – Semi-preparative HPLC of Methanolic Leaf Extract

ACE[®]
Ultra-inert
UHPLC & HPLC Columns

Application #AN3550

Conditions

Column:	ACE 5 C18
Dimensions:	250 x 7.75 mm
Part Number:	ACE-121-2508
Mobile Phase:	0.5% TFA in MeOH/H ₂ O (20:80 v/v)
Flow Rate:	4.5 mL/min
Injection:	500 µL
Temperature:	Ambient
Detection:	UV, 245 nm

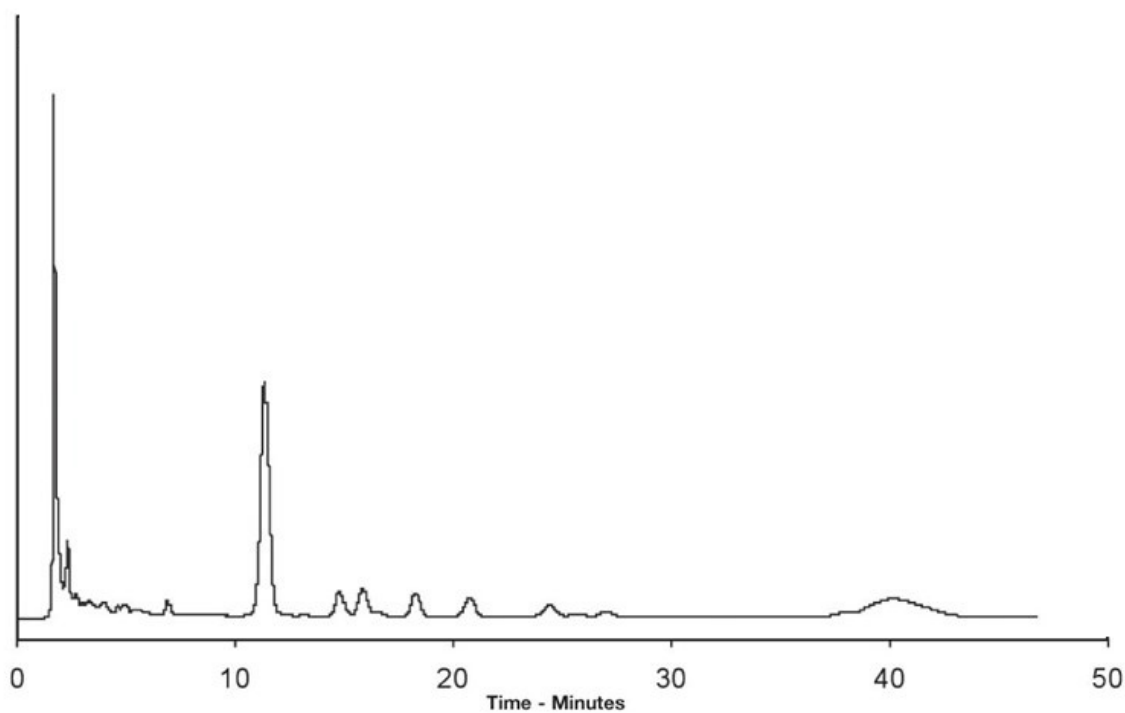

Reproduced with permission of School of Pharmacy, University of Sunderland, UK
© 2015 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com



ACE[®]
UHPLC & HPLC
Columns