



# Clopidogrel and Photodegradation Products

**ACE**<sup>®</sup>  
Ultra-inert  
UHPLC & HPLC Columns

Application #AN3110

## Conditions

Column:	ACE 5 C18
Dimensions:	150 x 4.6 mm
Part Number:	ACE-121-1546
Mobile Phase:	MeOH/aqueous TEA (pH 5.3 with H <sub>3</sub> PO <sub>4</sub> ) 75:25 v/v
Flow Rate:	1.2 mL/min
Injection:	20 µL
Temperature:	25 °C
Detection:	UV, 220 nm
Sample:	Exposed to UV light for 3.5 hours

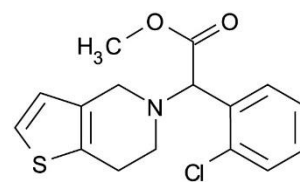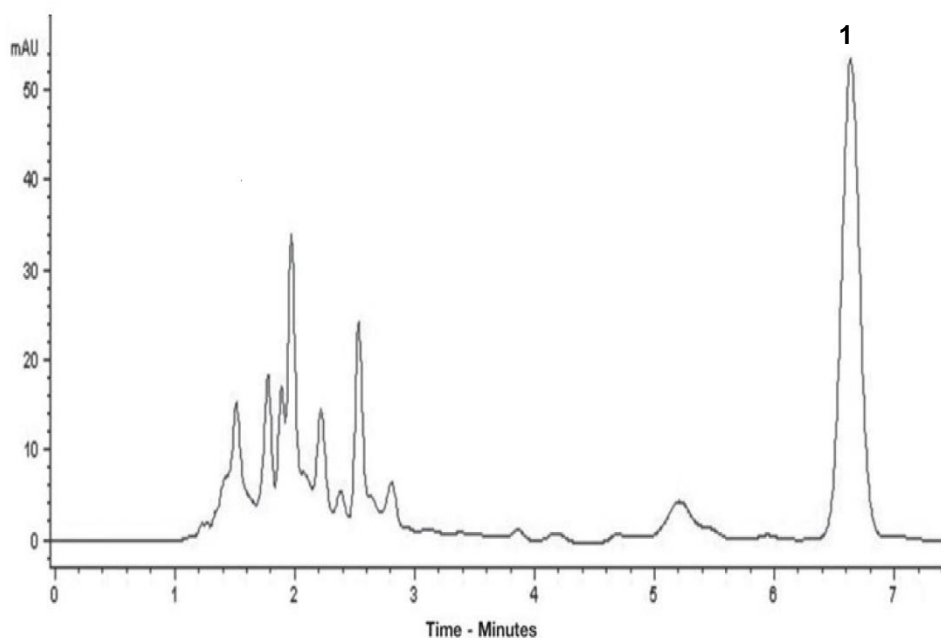

1. Clopidogrel

Reproduced with permission of Brazilian Pharmacopoeia  
© 2015 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit  
[www.ace-hplc.com](http://www.ace-hplc.com) or email: [info@ace-hplc.com](mailto:info@ace-hplc.com)

  
**ACE**<sup>®</sup>  
UHPLC & HPLC  
Columns