

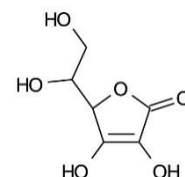


Conditions


Column:	ACE 5 C18
Dimensions:	250 x 4.6 mm
Part Number:	ACE-121-2546
Mobile Phase:	20 mM phosphoric acid in MeOH/H ₂ O (25:75 v/v)
Flow Rate:	1.5 mL/min
Injection:	5 µL
Temperature:	40 °C
Detection:	UV, 280 nm

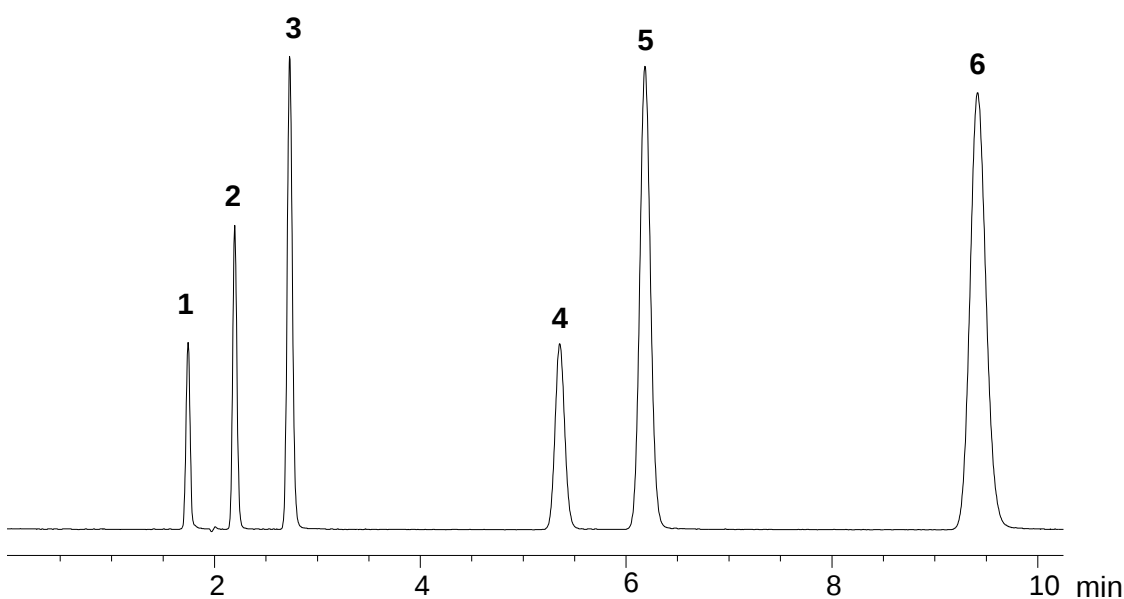


1. Ascorbic acid

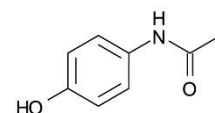


2. Paracetamol

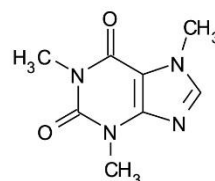


3. Caffeine

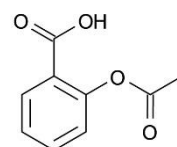

4. Aspirin

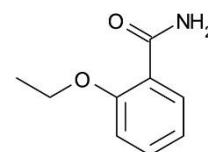

5. Ethenzamide

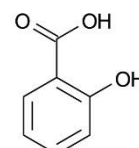

6. Salicylic acid