

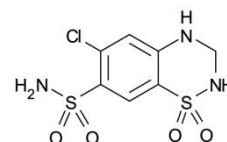


Conditions


Column:	ACE 5 CN
Dimensions:	200 x 4.6 mm
Part Number:	ACE-124-2046
Mobile Phase:	10 mM phosphoric acid in H ₂ O, pH 2.5/MeCN/MeOH (80:7:13 v/v/v)
Flow Rate:	1 mL/min
Injection:	20 µL
Temperature:	30 °C
Detection:	UV, 235 nm
Sample:	1 µg/mL each analyte

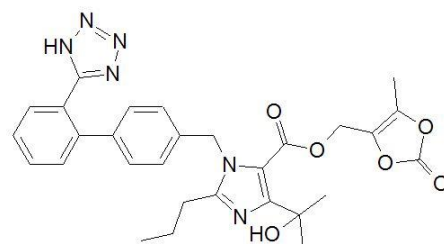


1. Hydrochlorothiazide

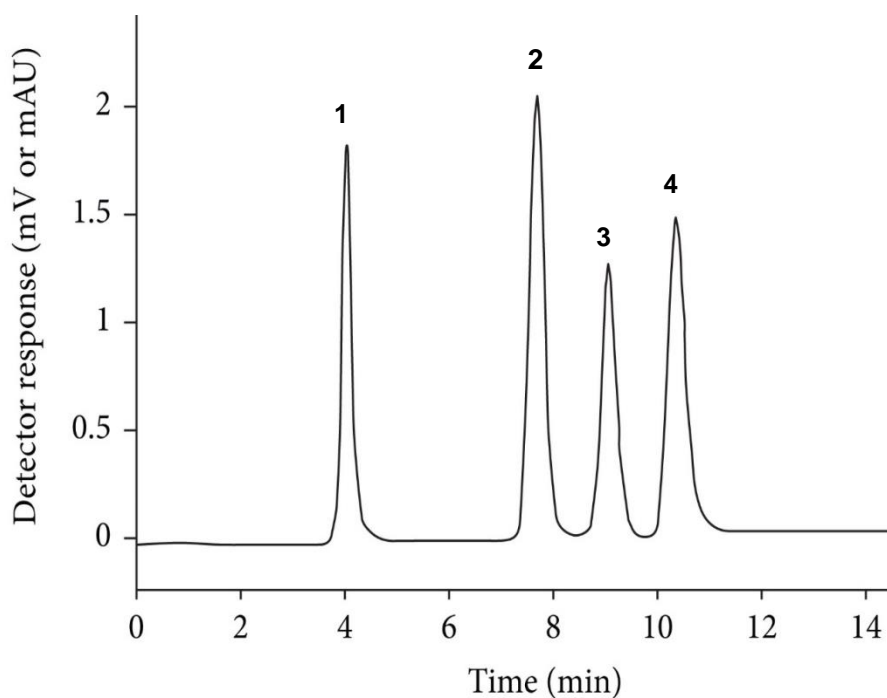

2. Olmesartan medoxomil

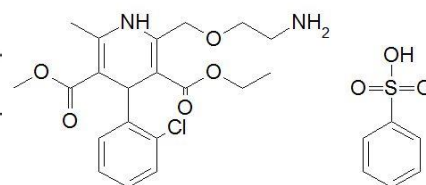

3. Amlodipine besylate

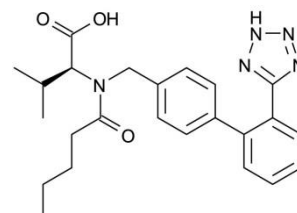

4. Valsartan

Tekkeli SEK. Development of an HPLC-UV Method for the Analysis of Drugs used for Combined Hypertension Therapy in Pharmaceutical Preparations and Human Plasma. Journal of Analytical Methods in Chemistry (2013)
<http://dx.doi.org/10.1155/2013/179627>

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.