


Determination of Amoxicillin Content

Application #AN6040

Conditions

Column: ACE Excel 5 SuperC18
Dimensions: 250 x 4.6 mm
Part Number: EXL-1211-2546U
Mobile Phase: 50 mM KH₂PO₄ (pH 5)/MeCN (97.5:2.5)
Flow Rate: 1.0 mL/min
Injection: 20 µL
Temperature: Ambient
Detection: UV, 254 nm
Monograph Method: ChP. 2015, Vol. 2

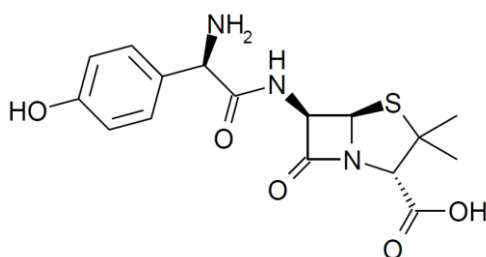

1. Amoxicillin

