


Determination of Amygdalin in Guizhi Fuling Capsule

ACE[®]
Ultra-inert
UHPLC & HPLC Columns

Application #AN5070

Conditions

Column: ACE Excel 5 C18
Dimensions: 250 x 4.6 mm
Part Number: EXL-121-2546U
Mobile Phase: MeOH/H₂O (20:80 v/v)
Flow Rate: 1.00 mL/min
Injection: 10 µL
Temperature: Ambient
Detection: UV, 218 nm
Method: ChP. 2015, Vol. 1

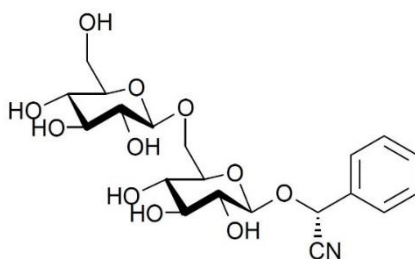

1. Amygdalin

