


Determination of Ephedrine Hydrochloride in Children's Chinese Medicine Prescription

ACE[®]
Ultra-inert
UHPLC & HPLC Columns

Application #AN5090

Conditions

Column: ACE Excel 5 SuperC18
Dimensions: 250 x 4.6 mm
Part Number: EXL-1211-2546U
Mobile Phase: MeCN/0.1 % phosphoric acid (aq.) (5:95 v/v)
Flow Rate: 1.00 mL/min
Injection: 10 µL
Temperature: Ambient
Detection: UV, 215 nm

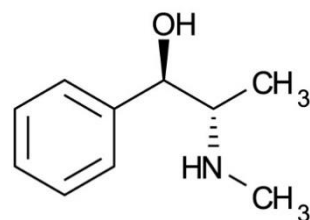

1. Ephedrine Hydrochloride

