


# Determination of Honokiol and Magnolol in Huoxiangzhengqi Water

**ACE**<sup>®</sup>  
Ultra-inert  
UHPLC & HPLC Columns


Application #AN6180

## Conditions

Column: ACE Excel 5 C18  
Dimensions: 250 x 4.6 mm  
Part Number: EXL-121-2546U  
Mobile phase: Methanol:Acetonitrile:Water (40:20:40 v/v/v)  
Flow Rate: 1.0 mL/min  
Injection: 10 µL  
Temperature: Ambient  
Detection: UV, 294 nm  
Monograph method: ChP. 2015, Vol. 1

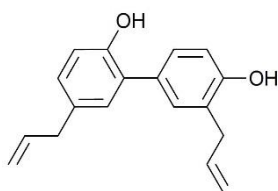

1. Honokiol

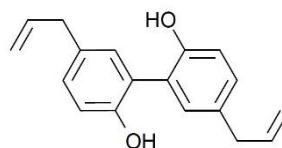

2. Magnolol

