
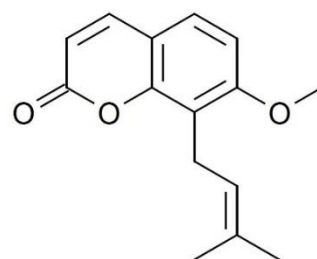



Determination of Isoimperatorin and Osthol in Tianma pills

Conditions

Column: ACE Excel 5 C18-PFP
Dimensions: 250 x 4.6 mm
Part Number: EXL-1210-2546U
Mobile Phase: MeCN/0.1 % phosphoric acid (aq.) (58:42 v/v)
Flow Rate: 1.0 mL/min
Injection: 10 µL
Temperature: 45 °C
Detection: UV, 320 nm
Monograph Method: ChP.2015, Vol. 1

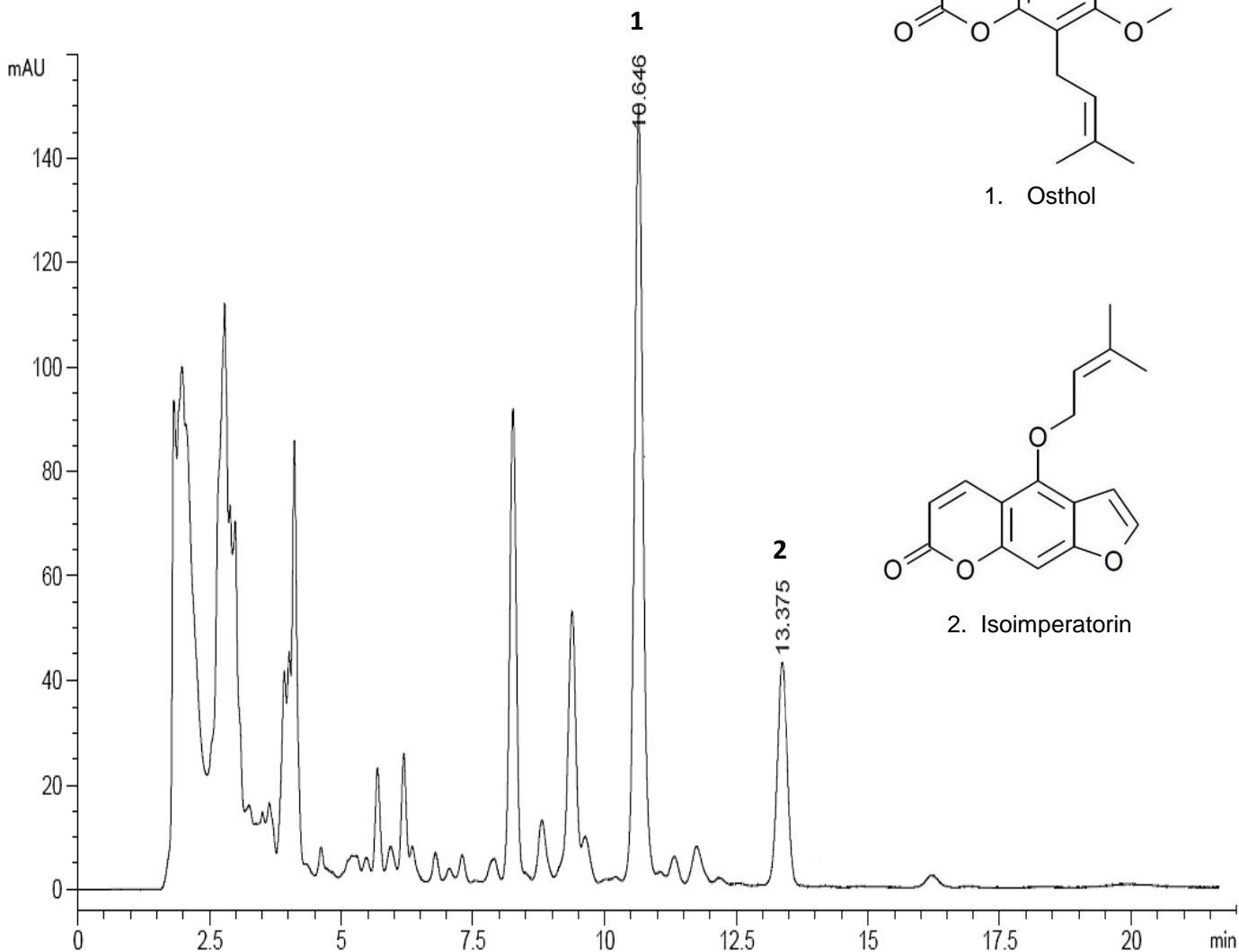

1. Osthol

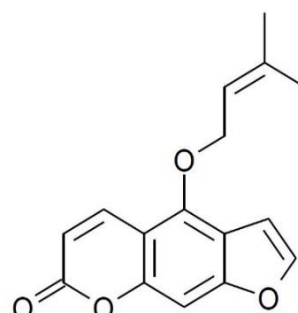

2. Isoimperatorin