


Determination of Loganin in Mingmu Dihuang Pills

ACE[®]
Ultra-inert
UHPLC & HPLC Columns


Application #AN6540

Conditions

Columns: ACE Excel 5 C18
Dimensions: 250 x 4.6 mm
Part Numbers: EXL-121-2546U
Mobile Phase: [THF/MeCN/MeOH (1:8:4 v/v/v)]/0.05 % H₃PO₄ solution (8:92 v/v)
Flow Rate: 1.0 mL/min
Injection: 10 µL
Temperature: Ambient
Detection: UV, 236 nm
Method: ChP.2015, Vol. 1

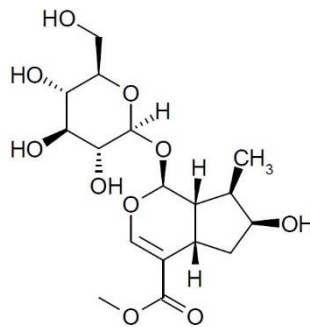

1. Loganin

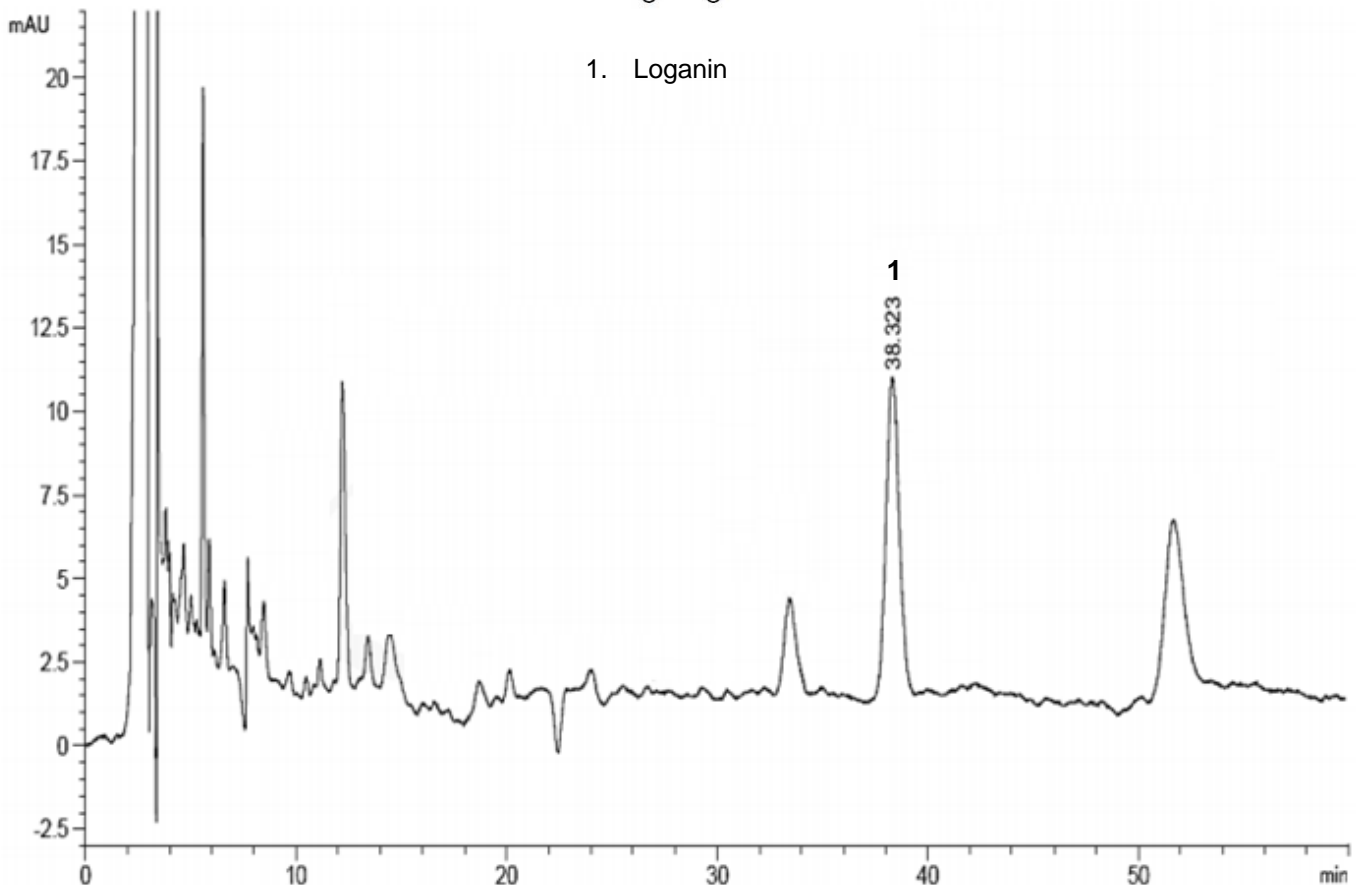

© 2018 ACT Ltd. All rights reserved. All products are available worldwide.


ACE[®]
UHPLC & HPLC
Columns

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com