


Determination of Paeoniflorin in Mingmu Dihuang Pills

ACE[®]
Ultra-inert
UHPLC & HPLC Columns

Application #AN6560

Conditions

Columns: ACE Excel 5 C18
Dimensions: 250 x 4.6 mm
Part Numbers: EXL-121-2546U
Mobile Phase: MeCN/0.1 % H₃PO₄ (aq) (14:86 v/v)
Flow Rate: 1.0 mL/min
Injection: 10 µL
Temperature: Ambient
Detection: UV, 230 nm
Method: ChP.2015, Vol. 1

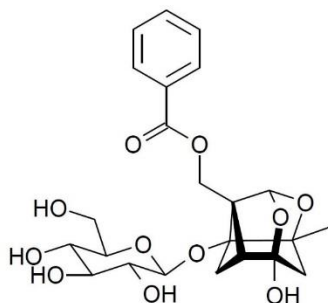

1. Paeoniflorin

