


# Determination of Paeoniflorin in Zhuifeng Tougu Capsules

**ACE**<sup>®</sup>  
Ultra-inert  
UHPLC & HPLC Columns


Application #AN6610

## Conditions

|               |                                   |
|---------------|-----------------------------------|
| Columns:      | ACE Excel 5 AQ                    |
| Dimensions:   | 150 x 4.6 mm                      |
| Part Numbers: | EXL-126-1546U                     |
| Mobile Phase: | MeOH/H <sub>2</sub> O (20:80 v/v) |
| Flow Rate:    | 1.0 mL/min                        |
| Injection:    | 5 µL                              |
| Temperature:  | Ambient                           |
| Detection:    | UV, 230 nm                        |

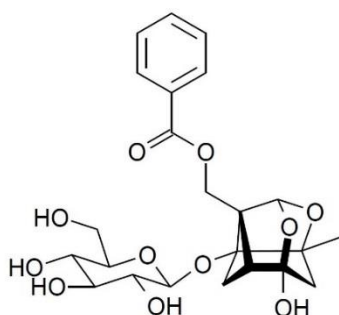

1. Paeoniflorin

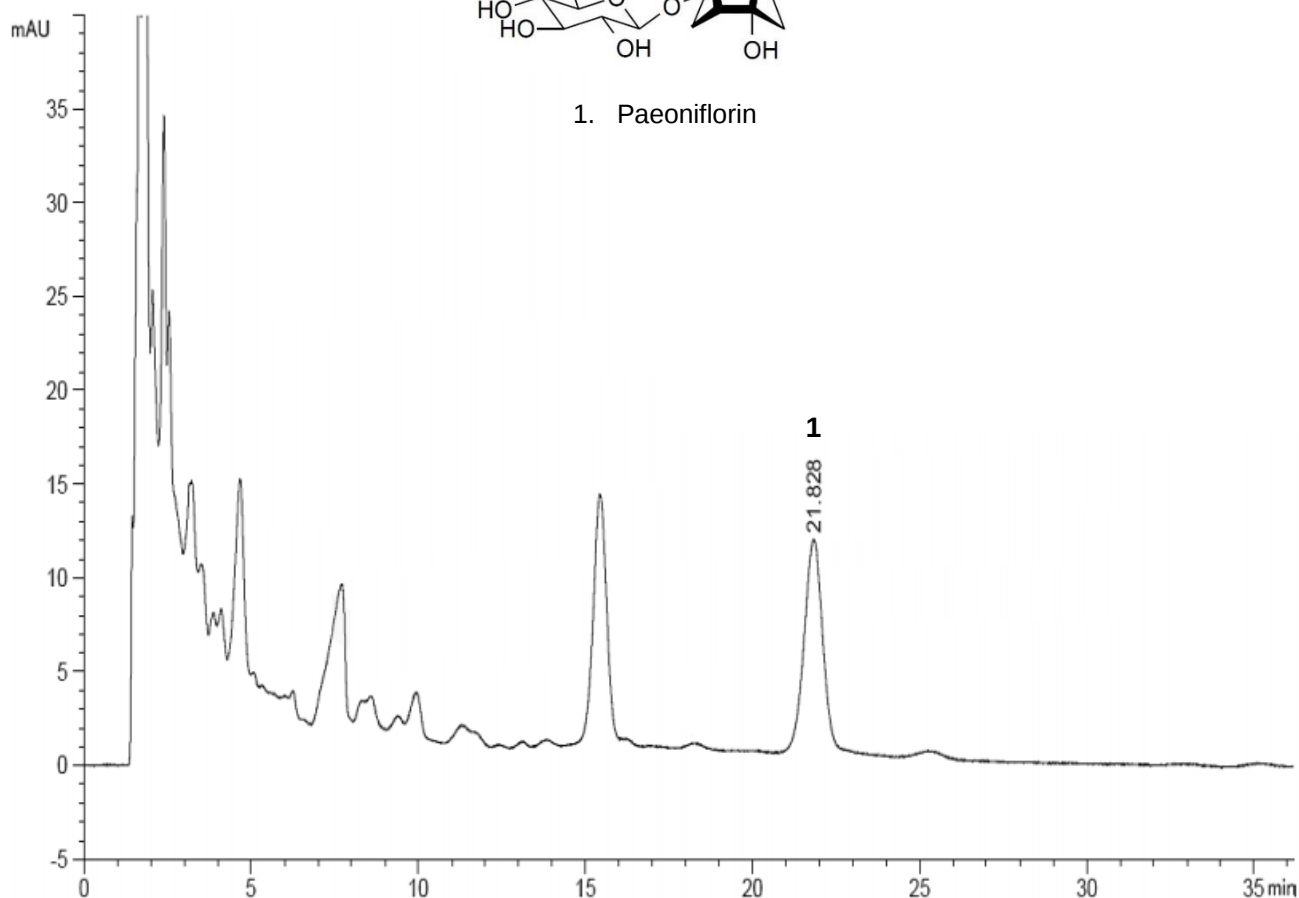

© 2018 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit  
[www.ace-hplc.com](http://www.ace-hplc.com) or email: [info@ace-hplc.com](mailto:info@ace-hplc.com)

  
**ACE**<sup>®</sup>  
UHPLC & HPLC  
Columns