
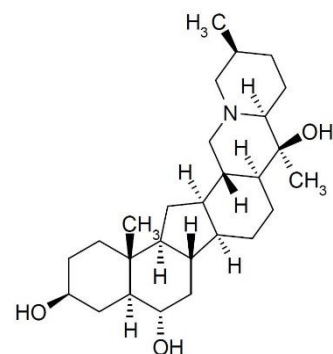


# Determination of Peimine and Peiminine in *Fritillaria thunbergii*


Application #AN6230

## Conditions

Column:	ACE Excel 5 SuperC18
Dimensions:	150 x 4.6 mm
Part Number:	EXL-1211-1546U
Mobile phase:	Acetonitrile:water:diethylamine (70:30:0.03 v/v/v)
Injection:	10 $\mu$ L
Temperature:	Ambient
Detection:	ELSD
Monograph method:	ChP. 2015, Vol. 1

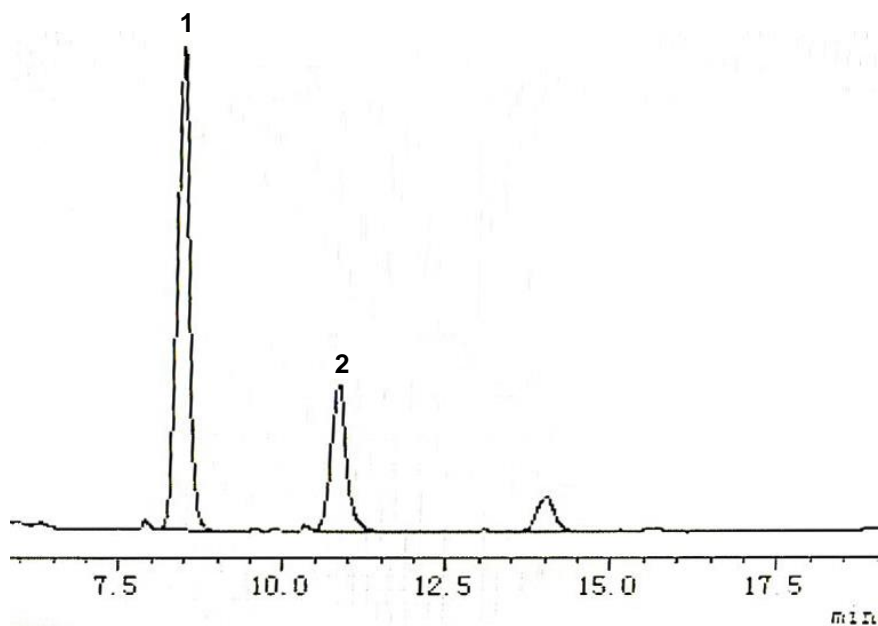

1. Peimine

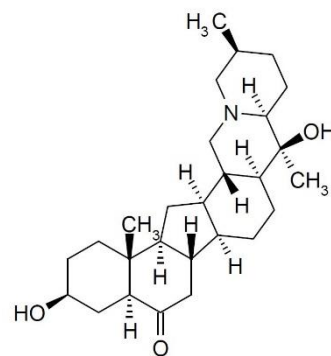

2. Peiminine