


Determination of the Content of Phenolic Ammonia Caffeine Tablets

ACE[®]
Ultra-inert
UHPLC & HPLC Columns


Application #AN6060

Conditions


Column: ACE Excel 5 C18
Dimensions: 250 x 4.6 mm
Part Number: EXL-121-2546U
Mobile Phase: 1 % acetic acid (pH 3.7)/MeOH (62:38 v/v)
Flow Rate: 1.0 mL/min
Injection: 10 µL
Temperature: Ambient
Detection: UV, 272 nm
Method: National Drug Standards (WS-10001-(HD-1021)-2002)

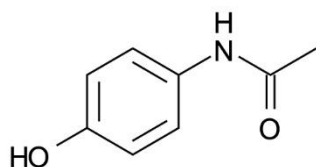


1. Acetaminophen (Paracetamol)

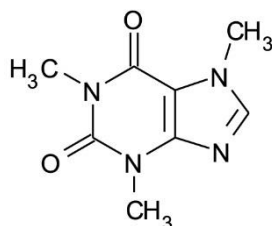

2. Caffeine

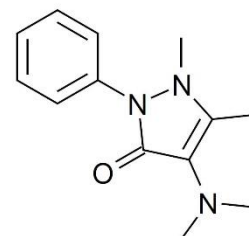

3. Aminopyrine

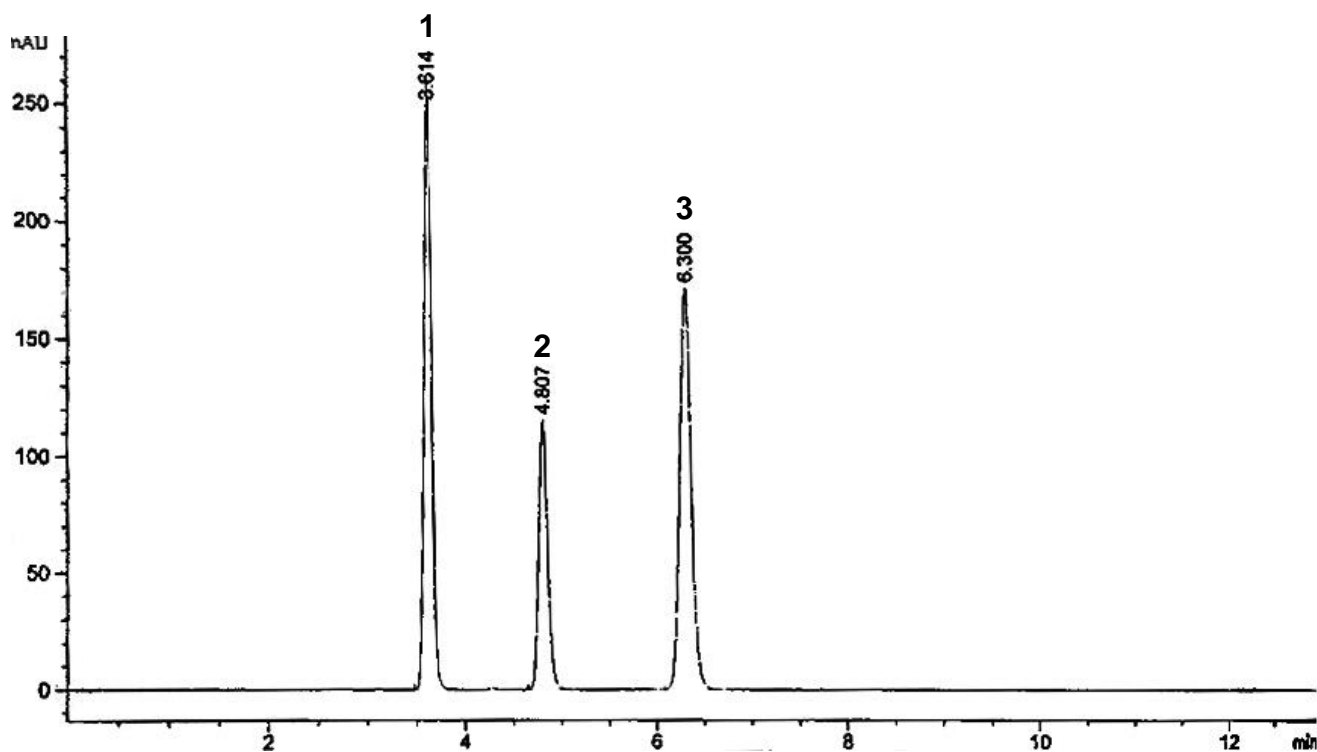

© 2018 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com


ACE[®]
UHPLC & HPLC
Columns