



Dexamethasone and Related Substances

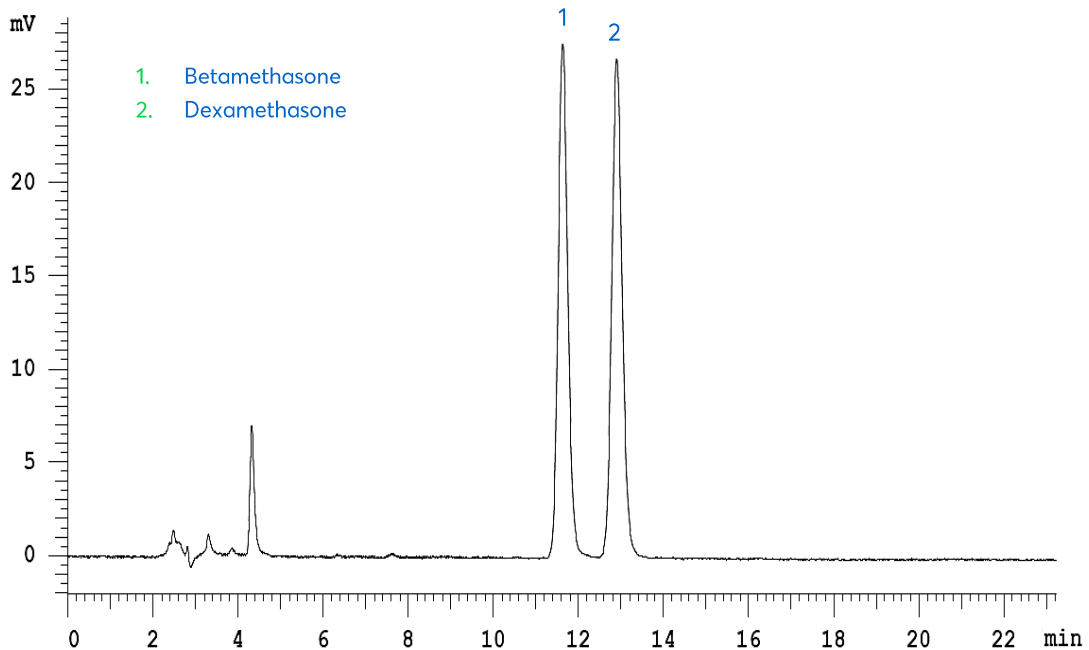


Method Details

CONDITIONS

Column:	ACE Excel 5 C18
Dimensions:	250 x 4.6 mm
Mobile Phase:	40 mM KH_2PO_4 and 24 mM hexylamine in $\text{MeCN}:\text{H}_2\text{O}$ (71:29 v/v)
Flow Rate:	1.0 mL/min
Injection:	20 μL
Temperature:	Ambient
Detection:	UV, 254 nm
Monograph method:	EP 6.0

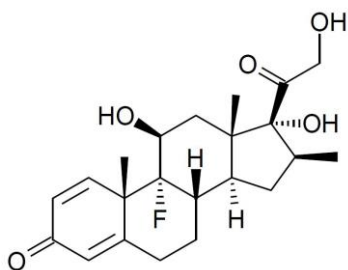

Structure 1: Betamethasone

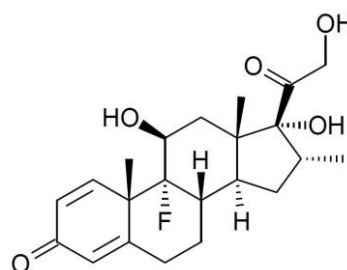

Structure 2: Dexamethasone

ORDERING TABLE

Product	Details	Size	Part Number
ACE Excel 5 C18	HPLC Column	250 x 4.6 mm	EXL-121-2546U