


EP Assay for Milk Thistle Dry Extract

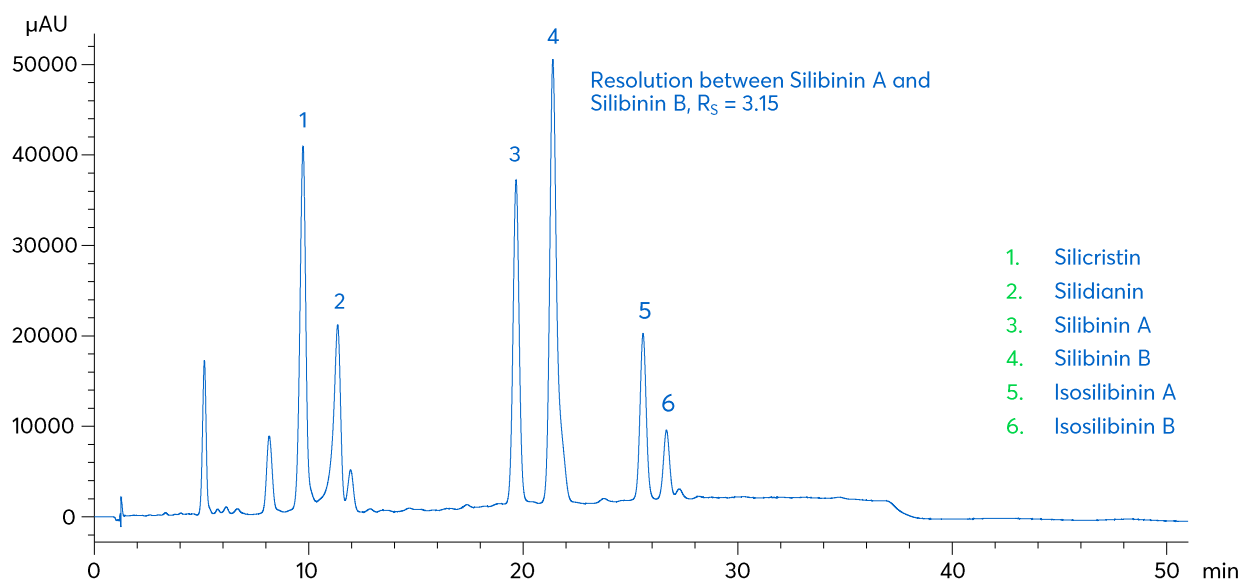

Method Details



CONDITIONS

Columns:	Avantor® ACE® 5 C18
Particle Size:	5 µm
Dimensions:	125 x 4.0 mm
Mobile Phase:	A: MeOH/H ₂ O/Orthophosphoric acid (35:65:0.5 v/v/v) B: MeOH/H ₂ O/Orthophosphoric acid (50:50:0.5 v/v/v)

Gradient:

Time (mins)	%B
0	0
28	100
35	100
36	0
51	0

Flow Rate:	0.8 mL/min
Injection Volume:	5 µL
Detection:	UV, 288 nm
Sample:	Milk thistle dry extract HRS
Method:	EP 10.0

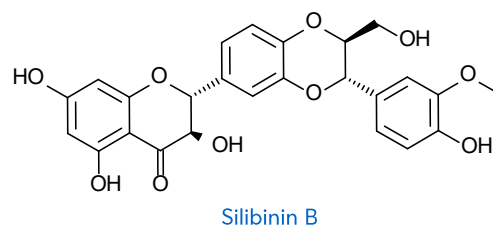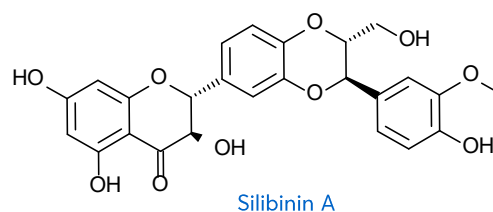

Always check the latest and official method information from the relevant pharmacopoeia prior to analysis

ORDERING TABLE

Product	Details	Size	Part number
Avantor® ACE® 5 C18	HPLC Column	125 x 4.0 mm	ACE-121-1204

Reproduced with permission of Stephan Avramov, PhD, Chromassist Ltd, Bulgaria