
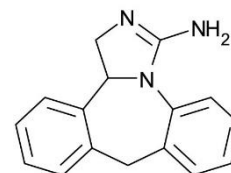


Conditions

Column:	ACE 5 C18
Dimensions:	150 x 4.6 mm
Part Number:	ACE-121-1546
Mobile Phase:	0.3% TEA pH 4.0 with phosphoric acid/MeOH (40:60 v/v)
Flow Rate:	1 mL/min
Injection:	20 µL
Temperature:	25 °C
Detection:	UV, 254 nm

(A) Epinastine Hydrochloride before UV radiation

(B) Epinastine Hydrochloride after 18 hours UV radiation

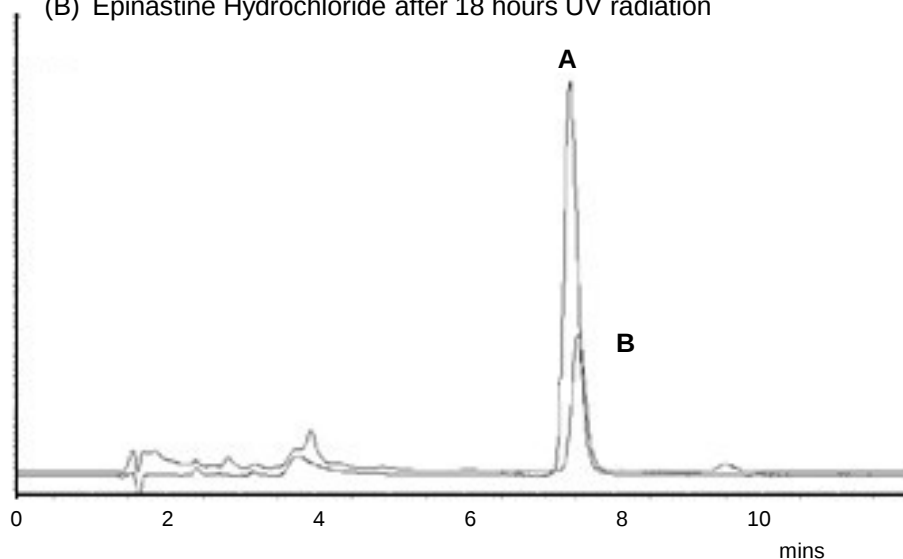

1. Epinastine

