

Ethanol Extract from Seed Cover (*Acacia Farnesiana*)

ACE[®]
Ultra-inert
UHPLC & HPLC Columns


Application #AN2900

Conditions

Column: ACE 5 C18
Dimensions: 250 x 4.6 mm
Part Number: ACE-121-2546
Mobile Phase: A: MeOH
B: H₂O

Time (mins)	%B
0.0	85
2.5	85
60.0	50
62.5	50
70.0	85

Flow Rate: 2 mL/min
Temperature: Ambient
Detection: UV, 230 nm

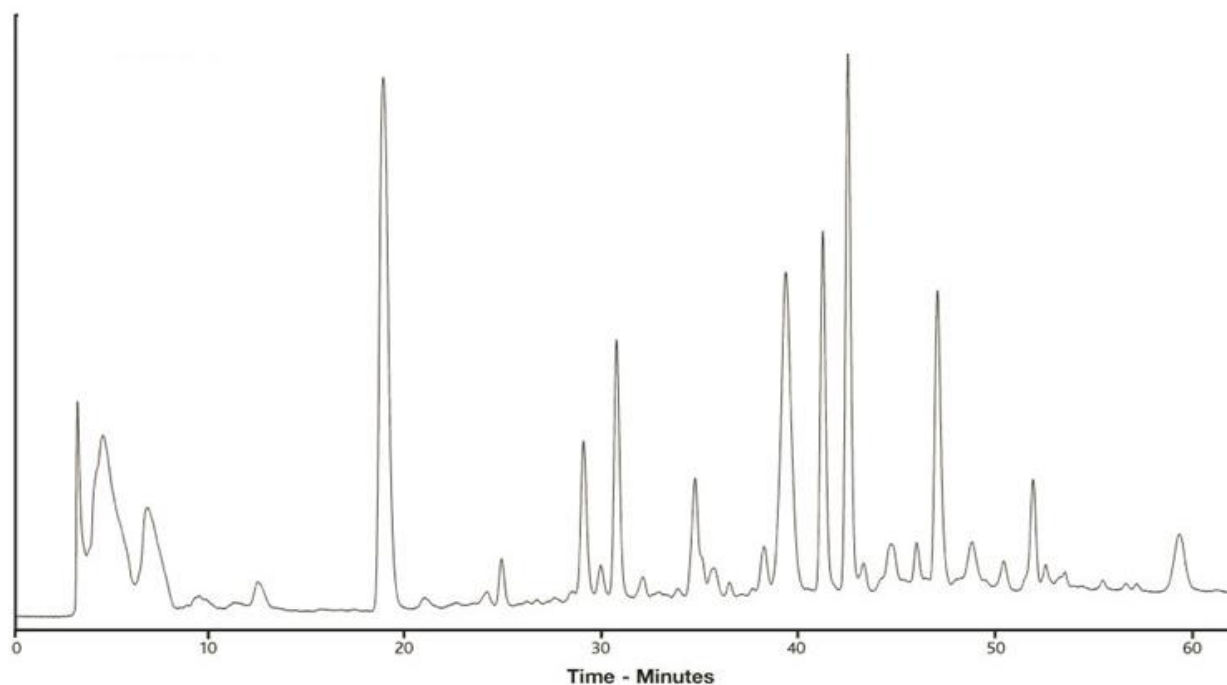

Reproduced with permission of School of Pharmacy, University of Sunderland, UK
© 2015 ACT Ltd. All rights reserved. All products are available worldwide.



For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com