


Garlic Analysis (I)

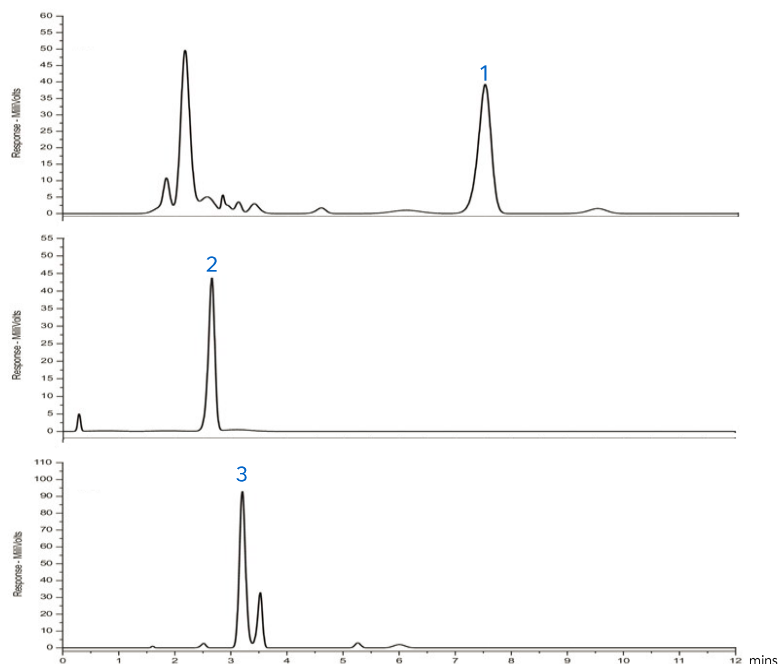


1. Allicin
2. Alliin
3. Deoxyalliin


Method Details

CONDITIONS


Column:	Avantor® ACE® 5 C18
Particle Size:	5µm
Dimensions:	250 x 4.6 mm
Mobile Phase:	H ₂ O/MeOH (50:50 v/v)
Flow Rate:	1 mL/min
Injection:	20 µL
Temperature:	30 °C
Detection:	UV, 210 nm

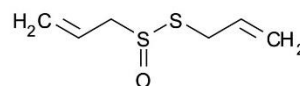

Allicin

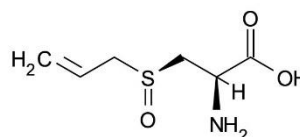

Alliin

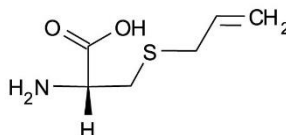

Deoxyalliin

ORDERING TABLE

Product	Details	Size	Part Number
Avantor® ACE® 5 C18	HPLC Column	250 x 4.6 mm	ACE-121-2546

Reproduced with permission of Neem Biotech, Cardiff, UK