
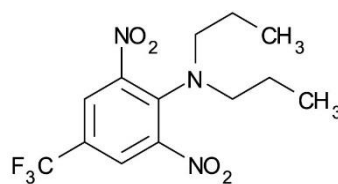


Herbicide - Trifluralin

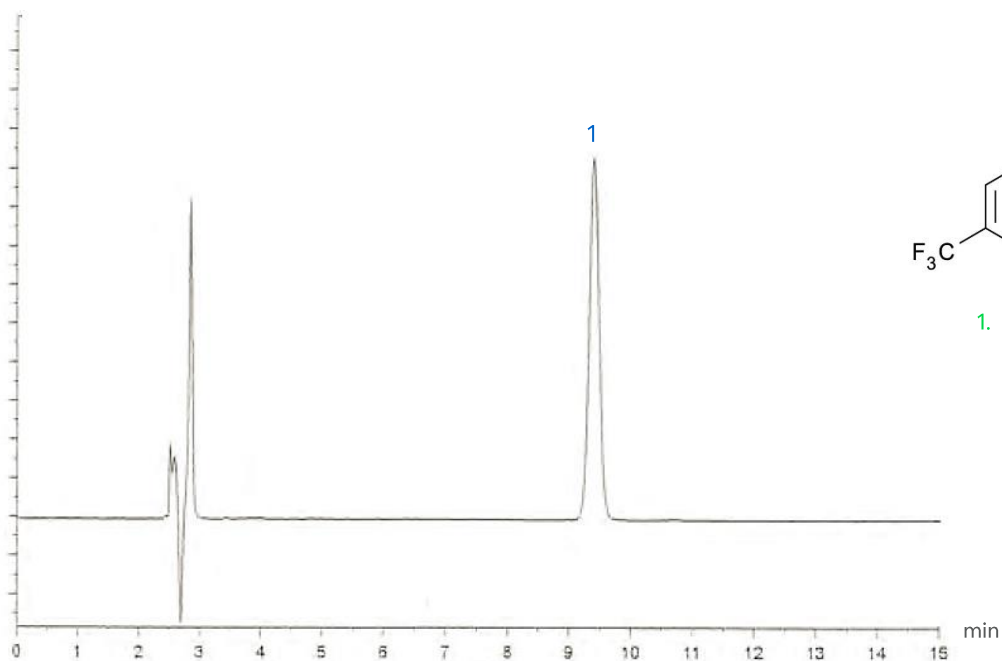

1. Trifluralin

Method Details

CONDITIONS

Column:	Avantor® ACE® 5 C18
Particle Size:	5 µm
Dimensions:	250 x 4.6 mm
Mobile Phase:	H ₂ O/MeOH (15:85 v/v)
Flow Rate:	1 mL/min
Temperature:	Ambient
Detection:	UV, 254 nm

ORDERING TABLE

Product	Details	Size	Part Number
Avantor® ACE® 5 C18	HPLC column	250 x 4.6 mm	ACE-121-2546
Methanol	VWR CHROMANORM® gradient grade	2.5 L	20864.320
Water	VWR CHROMANORM® super gradient grade	2.5 L	83650.320