


IP Method for Atorvastatin Calcium Assay


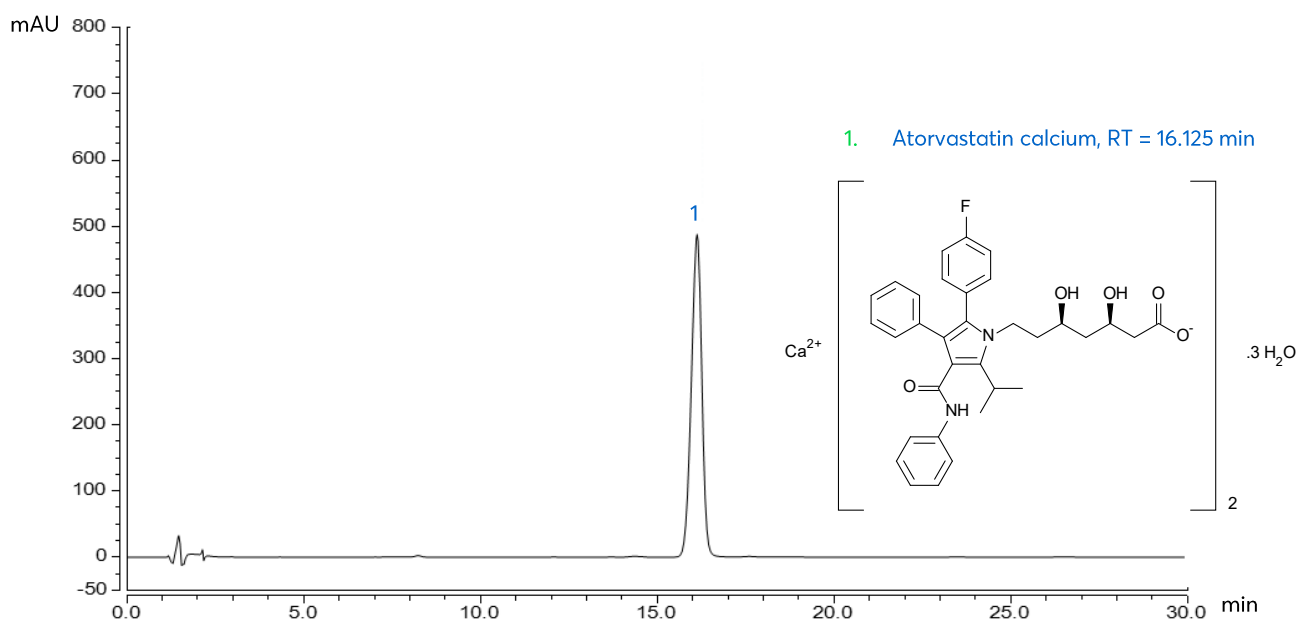


Method Details

CONDITIONS

Column:	Avantor® ACE® 5 C18
Particle Size:	5 µm
Dimensions:	250 x 4.6 mm
Solution A:	5.75 g ammonium dihydrogen orthophosphate in 1000 mL H ₂ O
Solution B:	MeCN/THF (92.5:7.5 v/v)
Mobile Phase:	Solution A/Solution B (58:42 v/v)
Flow Rate:	1.8 mL/min
Injection:	20 µL
Temperature:	Ambient
Detection:	UV, 246 nm
Method:	IP (2010)

Always check the latest and official method information from the relevant pharmacopoeia prior to analysis.

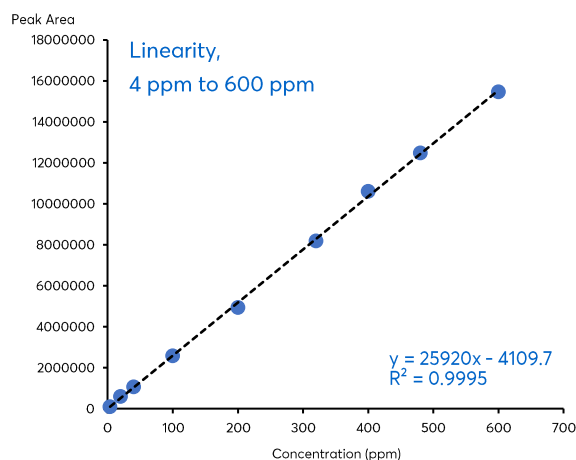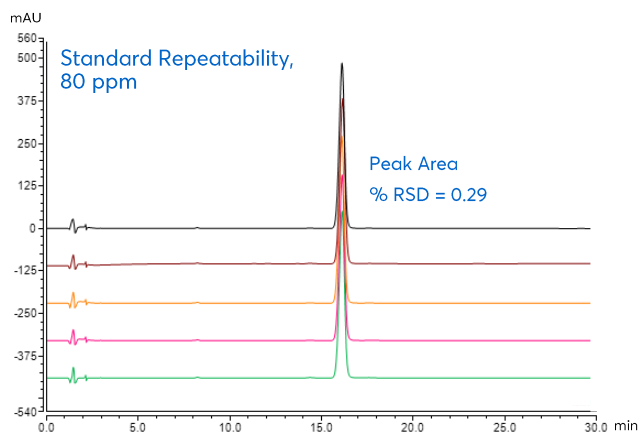

ORDERING TABLE

Product	Details	Size	Part number
Avantor® ACE® 5 C18	HPLC Column	250 x 4.6 mm	ACE-121-2546

Reproduced with permission of Mumbai Avantor Application lab.

Avantor® ACE®