



Ibuprofen in Combination with Antihistamine and Decongestant

ACE[®]
Ultra-inert
UHPLC & HPLC Columns


Application #AN2120

Conditions


Column: ACE Excel 3 C18-Amide
Dimensions: 150 x 4.6 mm
Part Number: EXL-1112-1546U
Mobile Phase: 0.01% potassium dihydrogen phosphate/MeCN (60:40 v/v)
Flow Rate: 0.6 mL/min
Injection: 0.5 µL
Temperature: 45 °C
Detection: UV, 214 nm

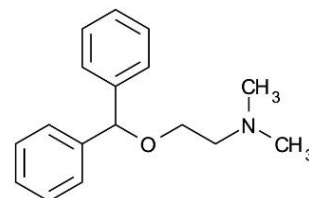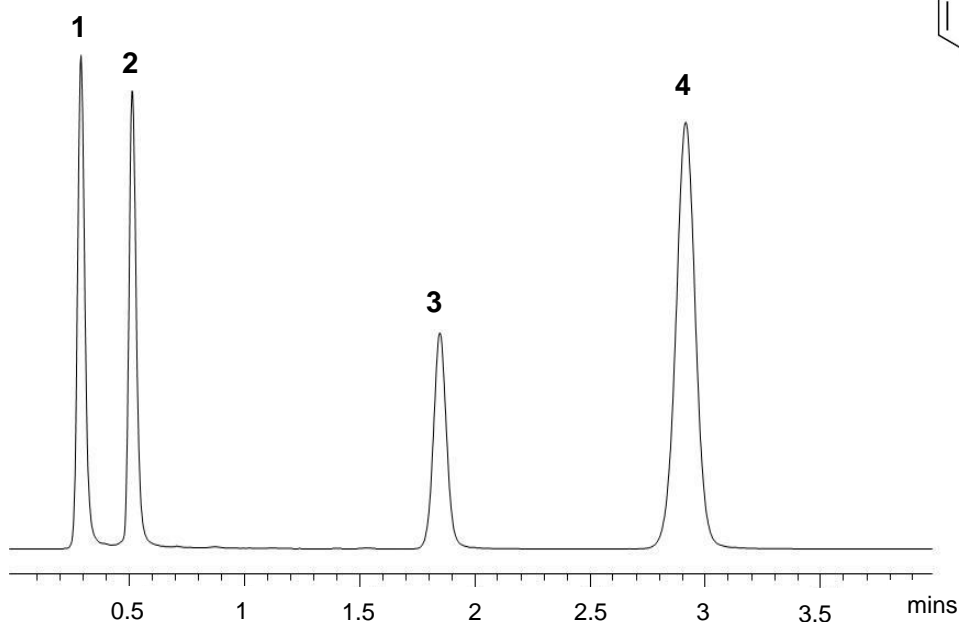


1. Diphenhydramine

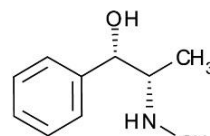

2. Pseudoephedrine

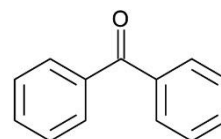

3. Benzophenone

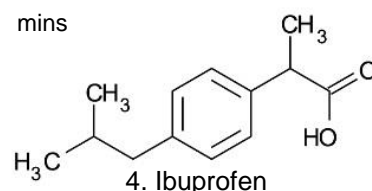

4. Ibuprofen