


### Conditions

Column: ACE 5 C18  
Dimensions: 250 x 4.6 mm  
Part Number: ACE-121-2546  
Mobile Phase: A: 0.1% TFA in H<sub>2</sub>O/MeCN (71:29 v/v)  
B: 0.1% TFA in H<sub>2</sub>O/MeCN (68:32 v/v)

| Time (mins) | %B |
|-------------|----|
| 0           | 10 |
| 16          | 90 |

Flow Rate: 1 mL/min  
Temperature: Ambient  
Detection: UV, 215 nm

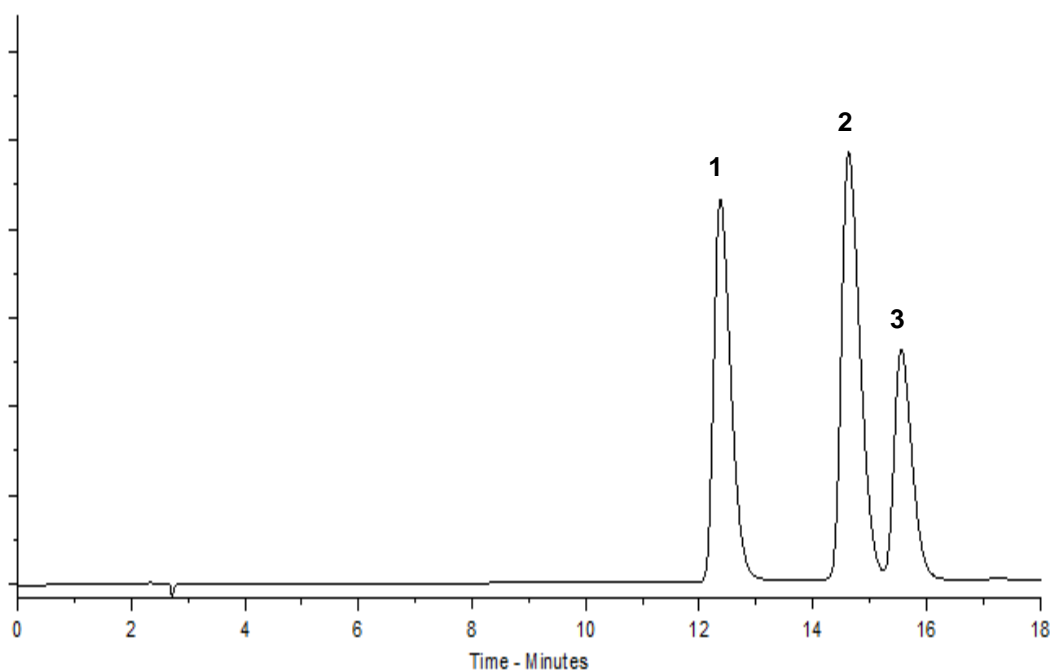

1. Bovine insulin
2. Human insulin
3. Porcine insulin