


Application note #7690

Analysis of Isavuconazonium

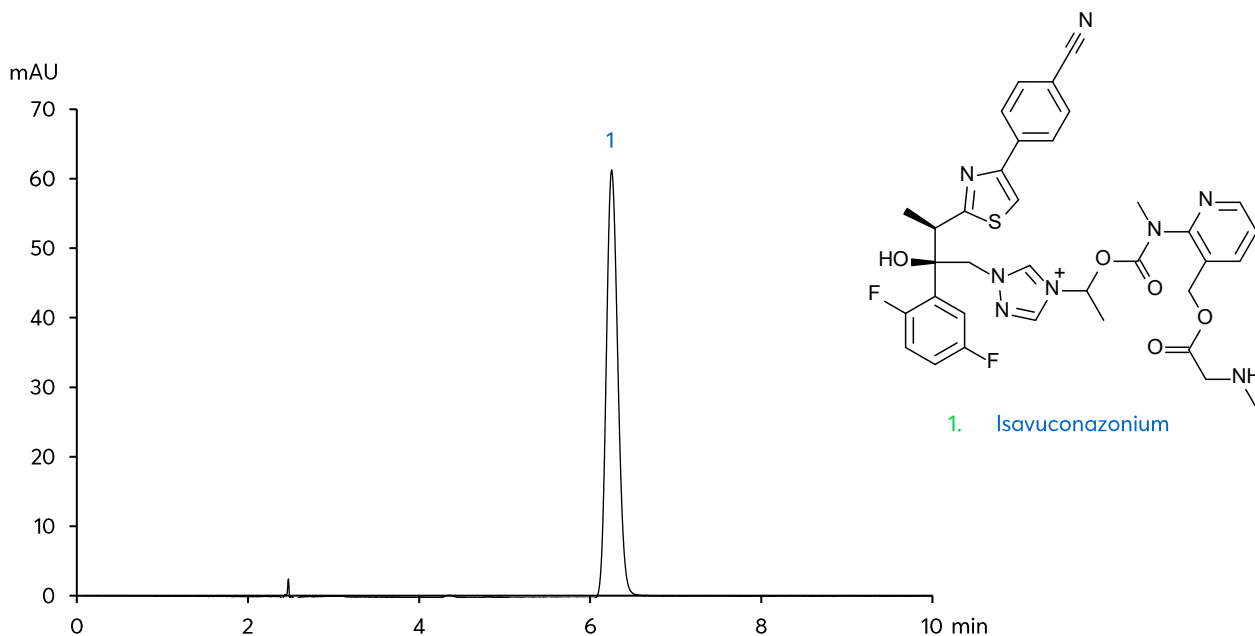

Method Details

CONDITIONS

Columns:	Avantor® ACE® Excel® 5 SuperC18
Particle Size:	5 µm
Dimensions:	250 x 4.6 mm
Mobile Phase:	0.1 M Ammonium acetate: MeOH (60:40)
Flow Rate:	1.0 mL/min
Injection Volume:	8 µL
Temperature:	30 °C
Detection:	UV, 280 nm

ORDERING TABLE

Product	Details	Size	Part Number
Avantor® ACE® Excel® 5 SuperC18	HPLC Column	250 x 4.6 mm	EXL-1211-2546U

Reproduced with permission of Mumbai Avantor Application lab.