

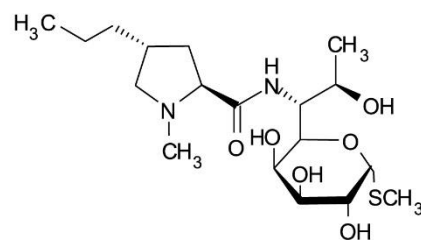


Conditions

Column:	ACE 5 C18
Dimensions:	250 x 4.6 mm
Part Number:	ACE-121-2546
Mobile Phase:	0.02 M sodium phosphate dibasic pH 3.0/MeCN (70:30 v/v)
Flow Rate:	1 mL/min
Injection:	25 µL
Temperature:	25 °C
Detection:	UV, 205 nm

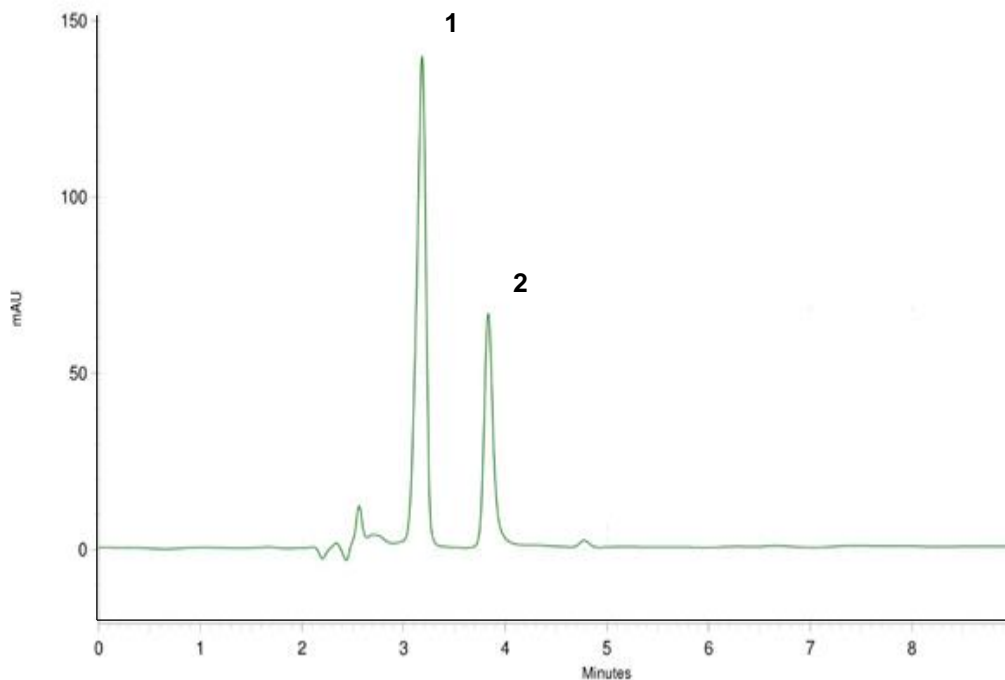

1. Lincomycin HCl

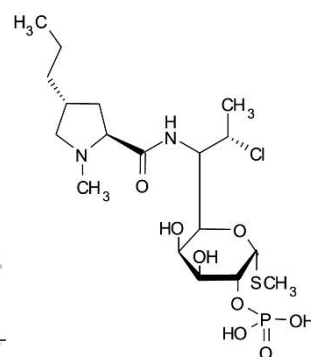

2. Clindamycin phosphate

Reproduced with permission of Department of Pharmacy, Faculty of Medicine and Surgery, University of Malta

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com