

Safety Funnel with Hinged Lid

Strong Protection when collecting
Liquid Waste!



Clever Technology, Strong Protection!

Safety Funnel

✓ Improved Closing

Comfortable handling: easy, quick opening and closing with "click"-closure.

✓ Hinged Lid with Sealing Ring

For secure closure, when the funnel is not in use.

✓ Stable Hinge

A double-screw hinge, constructed of rust-free stainless steel: for long life and stability.

✓ Optimized Dirt Sieve

For collecting mixing magnets or large dirt particles. Laboratory bottles can also be placed here to dry off. The hole size has been optimized, in order to improve the running-off of chemicals. The sieve can be removed for cleaning.

✓ Various Screw-Threads



Suitable for different sizes of collecting canister (see Table on Page 3).

✓ Safety Lance



Ensures for safe and good flowing of sensitive chemicals, without any "stuttering."


Dimensions

*Installation height:

+4.13 in

Lid open vertical:

+10.63 in

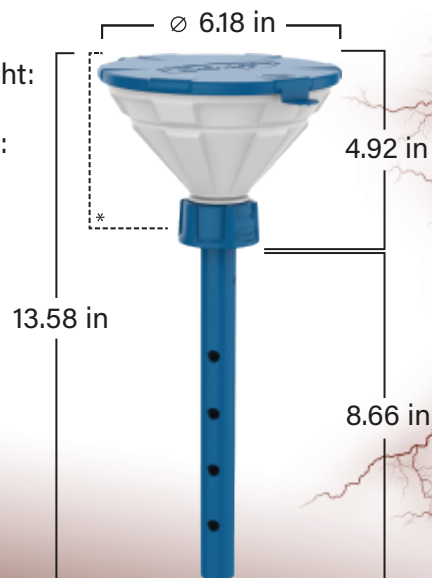

Configurations



Fig.	Part No.	Screw-Thread Of Container	Lance	Electrically conductive (incl. grounding cable and clamp)	Integrated Level Control (with Float in mm)	Delivery with Adapter
A	MM318 952	GL 45	-	-	-	-
-	MM318 962	GL 45	-	●	-	-
B	MM318 992	GL 45	●	-	-	-
-	MM318 955	S 50	-	-	-	●
-	MM318 985	S 50	●	●	-	●
-	MM318 995	S 50	●	-	-	●
-	MM318 953	S 51	-	-	-	●
-	MM318 983	S 51	●	●	-	●
-	MM318 993	S 51	●	-	-	●
-	MM318 951	S 55	-	-	-	●
-	MM318 961	S 55	-	●	-	●
-	MM318 981	S 55	●	●	-	●
-	MM318 991	S 55	●	-	-	●
-	MM318 950	S 60/61	-	-	-	-
C	MM318 960	S 60/61	-	●	-	-
D	MM318 980	S 60/61	●	●	-	-
-	MM318 990	S 60/61	●	-	-	-
E	MM317 637	S 60/61	-	●	● 150	-
-	MM317 651	S 60/61	-	●	● 190	-
-	MM318 954	S 65	-	-	-	●
-	MM318 964	S 65	-	●	-	●
-	MM318 984	S 65	●	●	-	●
-	MM318 994	S 65	●	-	-	●



Information about "Delivery with Adapter:"

The safety funnel is available in various container screw-thread sizes. In order to be environmentally-friendly, certain variants are delivered together with a screw-thread adapter. Those to which this applies can be seen in the column "Delivery with Adapter."





**THE No. 1 System
for Solvent Safety**

MAC-MOD Analytical

103 Commons Court, Chadds Ford PA 19317



610-358-9696



efaden@mac-mod.com



mac-mod.com