

# Non-Steroidal Anti-Inflammatory Drugs (III)


Application #AN3570

## Conditions

Column: ACE UltraCore 2.5 SuperC18  
Dimensions: 50 x 3.0 mm  
Part Number: CORE-25A-0503U  
Mobile Phase: A: 0.1% formic acid in H<sub>2</sub>O  
B: 0.1% formic acid in MeCN

Time (mins)	%B
0.02	20
2.71	70
3.39	70
3.52	20

Flow Rate: 0.85 mL/min  
Injection: 1.04 µL  
Temperature: 40 °C  
Detection: UV, 254 nm

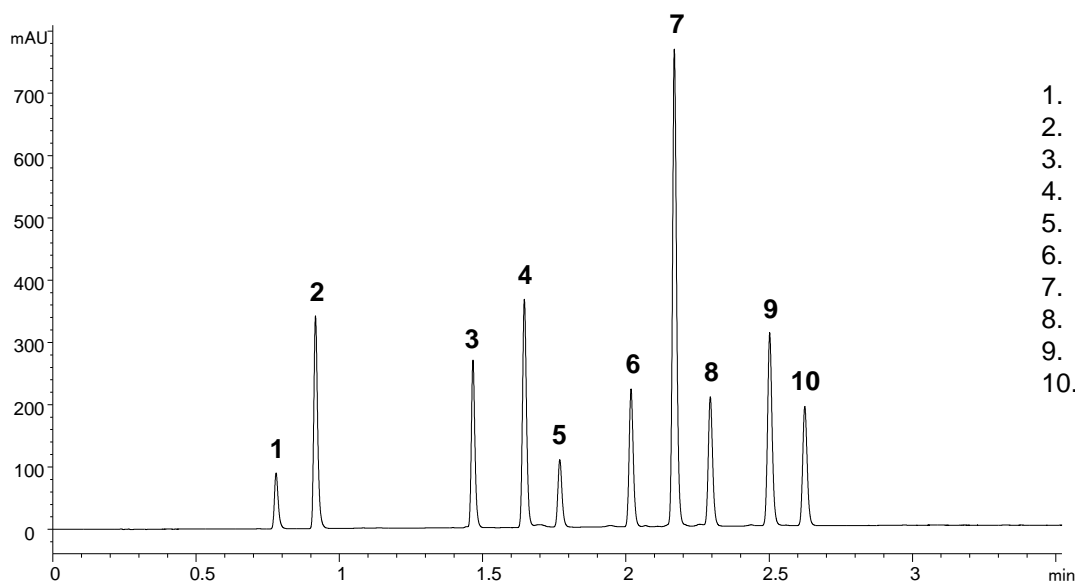

1. Aspirin
2. Phenacetin
3. Sulindac
4. Tolmetin
5. Naproxen
6. Nimesulide
7. Flurbiprofen
8. Diclofenac
9. Phenylbutazone
10. Meclofenamic acid