

Notoginseng Triterpenes in Xueshuantong Injection on ACE UltraCore 2.5 SuperC18

ACE[®]
Ultra-inert
UHPLC & HPLC Columns

Application #AN6500


Conditions

Column: ACE UltraCore 2.5 SuperC18
Dimensions: 100 x 4.6 mm
Part Number: CORE-25A-1046U
Mobile Phase: A: H₂O
B: MeCN

Time (mins)	%B
0	18
20	18
45	46
46	95
48	95
49	18
70	18

Flow Rate: 1.0 mL/min
Injection: 10 µL
Temperature: 30 °C
Detection: UV, 203 nm

1. Notoginsenoside R1
2. Ginsenoside Rg1
3. Ginsenoside Re
4. Ginsenoside Rb1
5. Ginsenoside Rd

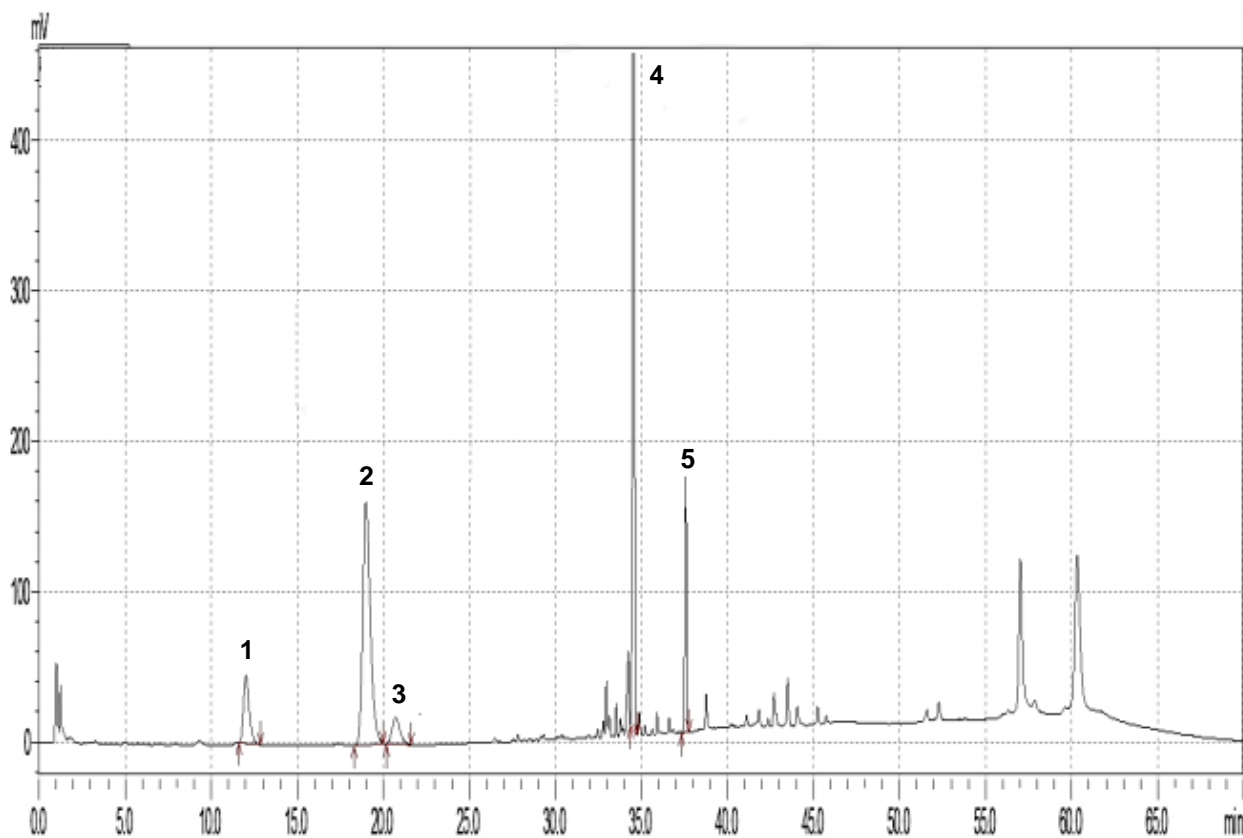

© 2018 ACT Ltd. All rights reserved. All products are available worldwide.



For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com