


Pantoprazole – (Protium)

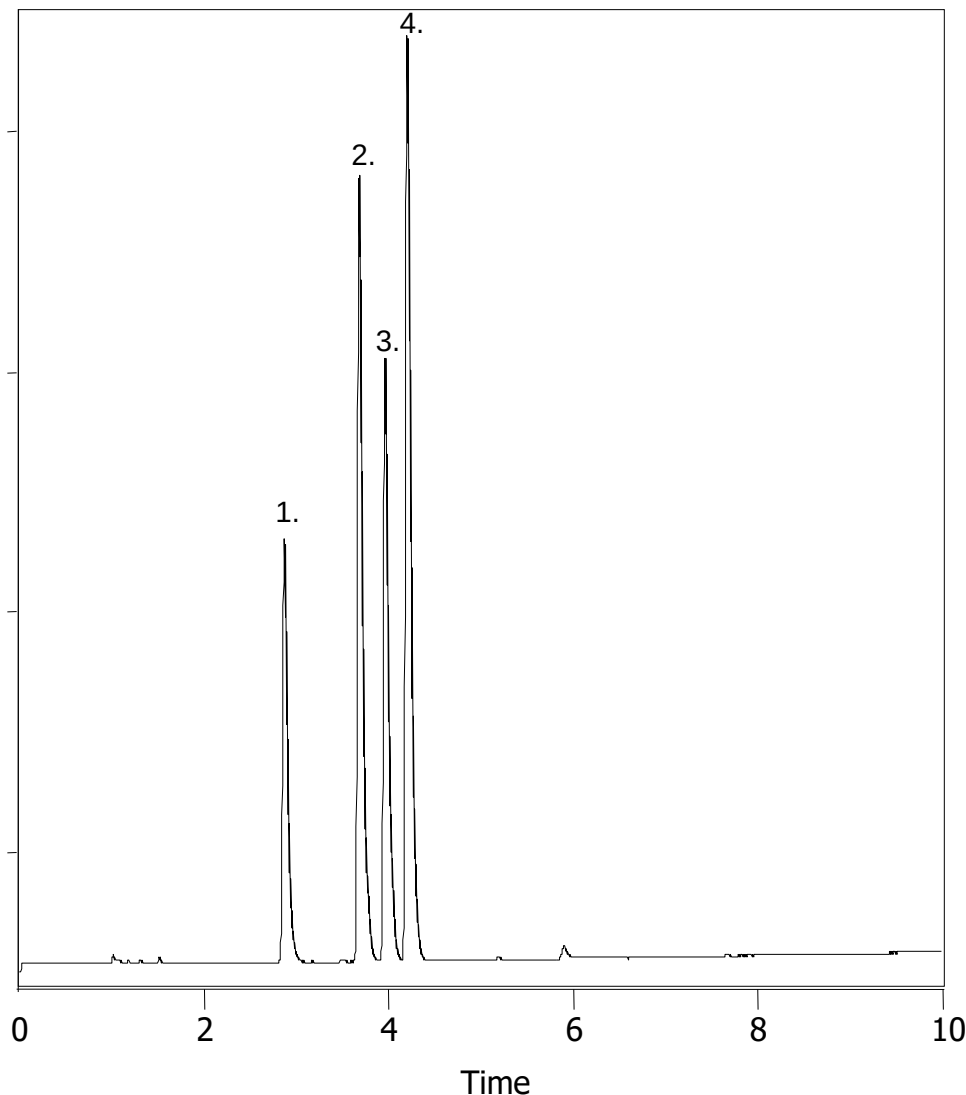

Column : **2.6μ SpeedCore RP Amide
150x4.6mm**

p/n : **SPRA-050726**

Mobile Phase:

A: H₂O + 0.1% Formic acid

B: ACN + 0.1% Formic acid

Gradient: **20-50% B in 10 minutes**

Flow : **1.0ml/min**

Temp : **30°C**

λ : **254nm**

- 1. Dopamine**
- 2. Imipramine**
- 3. Nortriptyline**
- 4. Trimipramine**