

Parotoid Macroglad Secretions from South American Toads

ACE[®]
Ultra-inert
UHPLC & HPLC Columns


Application #AN3970

Conditions

Column: ACE 5 C18
Dimensions: 250 x 4.6 mm
Part Number: ACE-121-2546
Mobile Phase: A: 0.1% TFA in H₂O
B: 0.1% TFA in MeCN/H₂O (90:10 v/v)


Time (mins)	%B
0	0
5	0
25	100

Flow Rate: 1 mL/min
Detection: PDA, 214 nm (Scanning 200-500 nm)

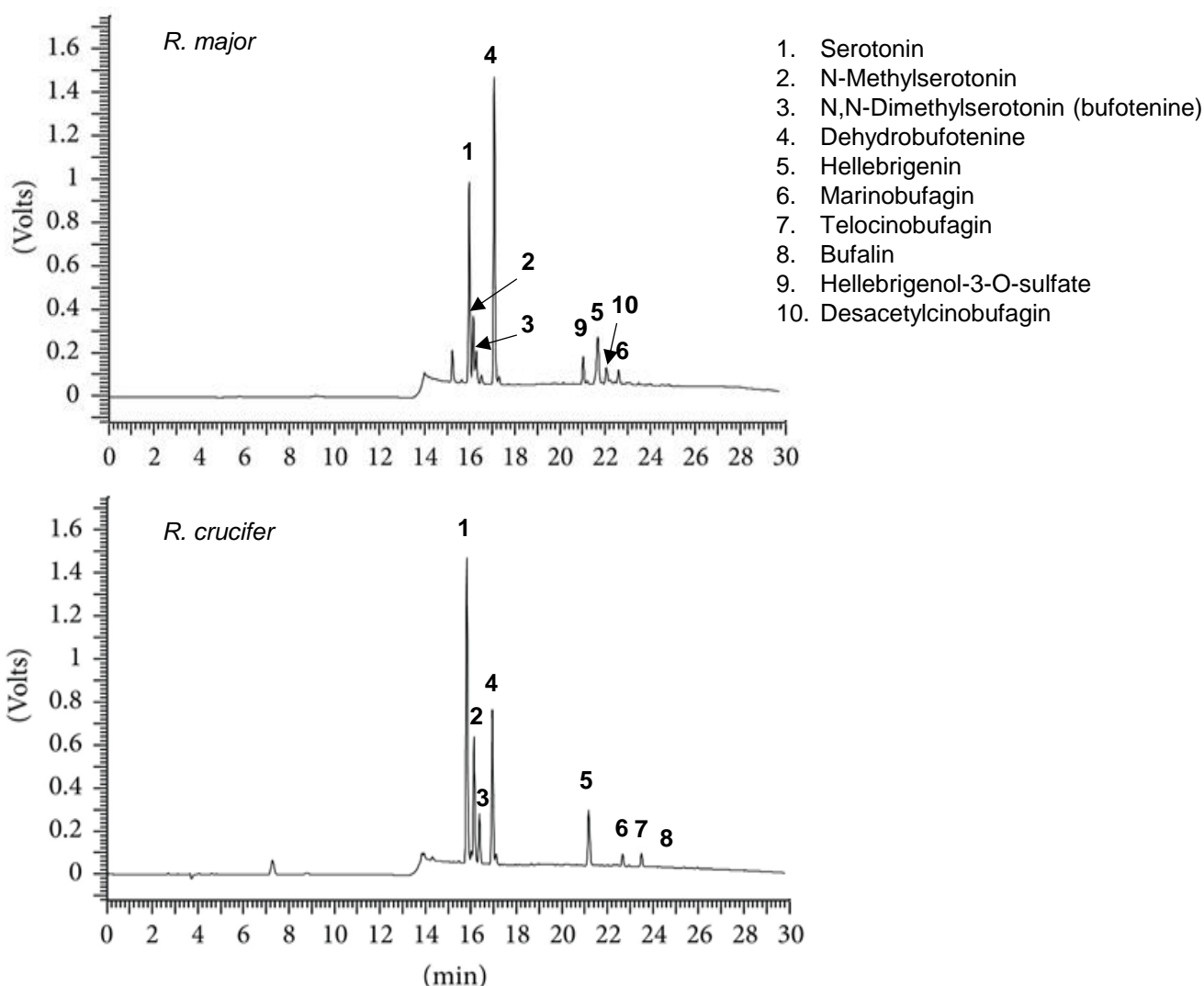

Sciani JM, Angeli CB, Antoniazzi MM, Jared C, Pimenta DC. Differences and Similarities among Parotoid Macroglad Secretions in South American Toads: A Preliminary Biochemical Delineation. The Scientific World Journal 2013, <http://dx.doi.org/10.1155/2013/937407>

© 2015 ACT Ltd. All rights reserved. All products are available worldwide.

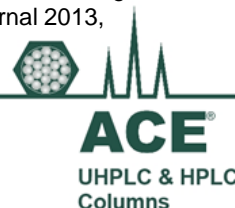

For more information contact your local ACE distributor or visit
www.ace-hplc.com or email: info@ace-hplc.com