


Phenones



- Column : **Evosphere AQUA 150x4.6 3um**
- p/n : **EVOAQ-050703**
- Mobile Phase:
 - **A- 0.1% Formic Acid in Water**
 - **B- 0.1% Formic Acid in MeOH**
 - **0-40% B in 15 mins**
- Flow : **1.0 ml/min**
- Temp : **20°C**
- λ : **254nm**

1. **Acetophenone**
2. **Propiophenone**
3. **Butyrophenone**
4. **Benzophenone**
5. **Valerophenone**

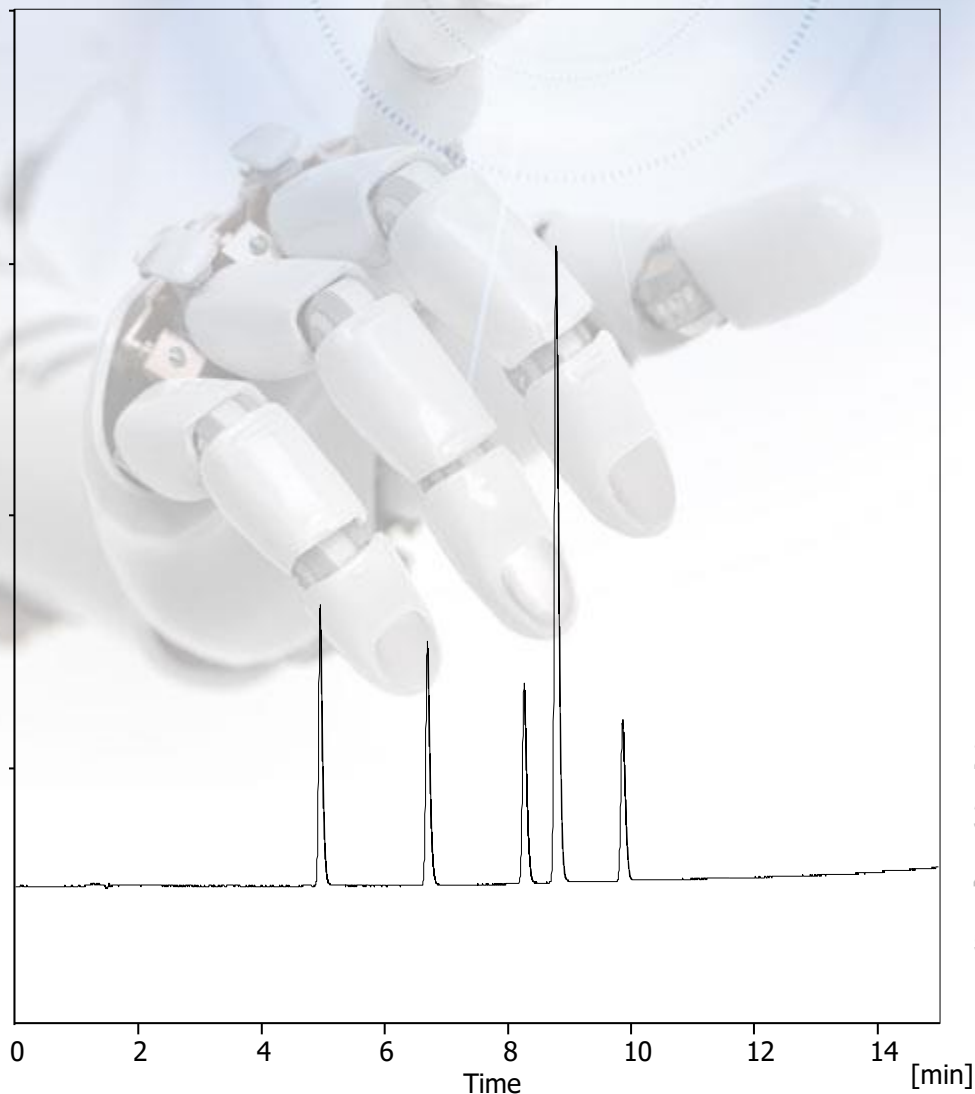