


# Phenoxy Acids

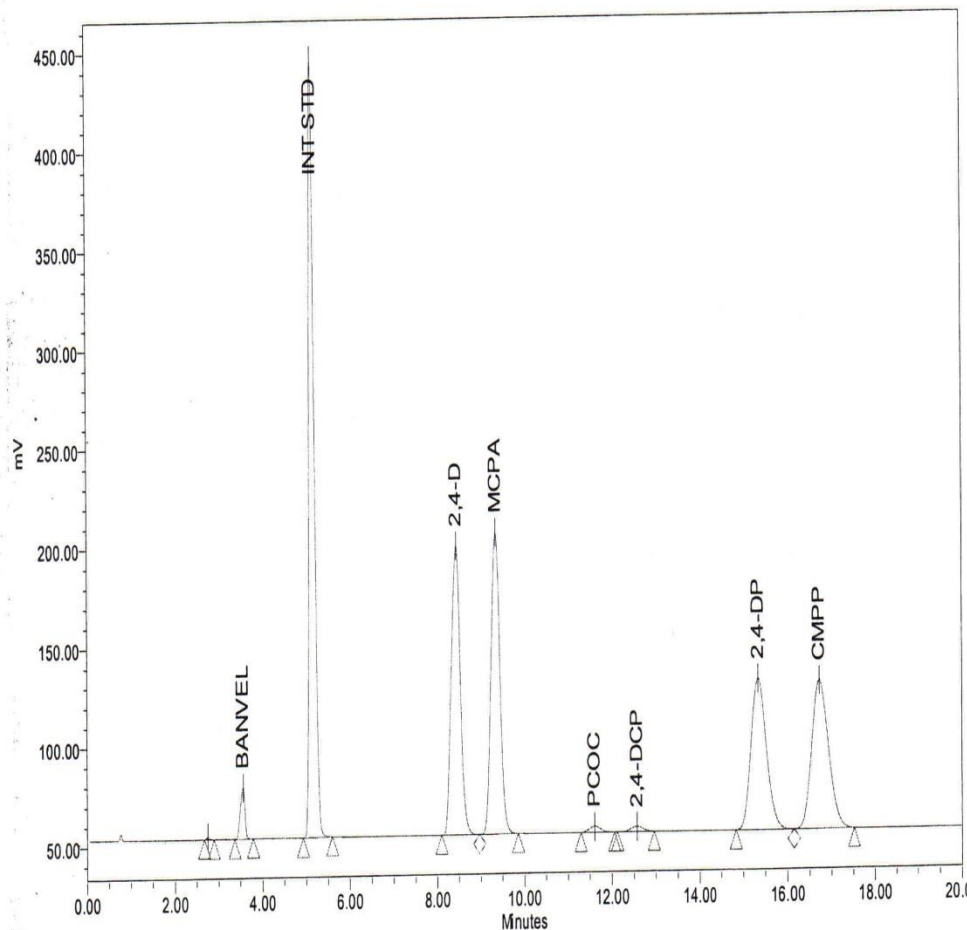

Column : **3 $\mu$  Fortis C18**

**100x4.6mm**

p/n : **F18-050503**

Mobile Phase:

**70:30 Water Ac(0.1%):ACN**

Flow : **1.5ml/min**

Temp : **50°C**

$\lambda$  : **254nm**

- 1. Banvel**
- 2. Internal STD**
- 3. 2,4 D**
- 4. MCPA**
- 5. PCOC**
- 6. 2,4-DCP**
- 7. 2,4-DP**
- 8. CMPP**