


Phthalate Selectivity



- Column : **Evosphere AQUA 150x4.6 3um**
- p/n : **EVOAQ-050703**
- Mobile Phase:
 - A- Water**
 - B- ACN**
 - 10-100% B in 20 mins, 10 min hold**
- Flow : **1.0 ml/min**
- Temp : **AMB**
- λ : **210nm**

- 1. Benzyl butyl phthalate**
- 2. Bis(2-ethylhexyl) adipate**
- 3. Bis(2-ethylhexyl) phthalate**
- 4. Dibutyl phthalate**
- 5. Diethyl phthalate**
- 6. Dimethyl phthalate**
- 7. Di-*n*-octyl phthalate**

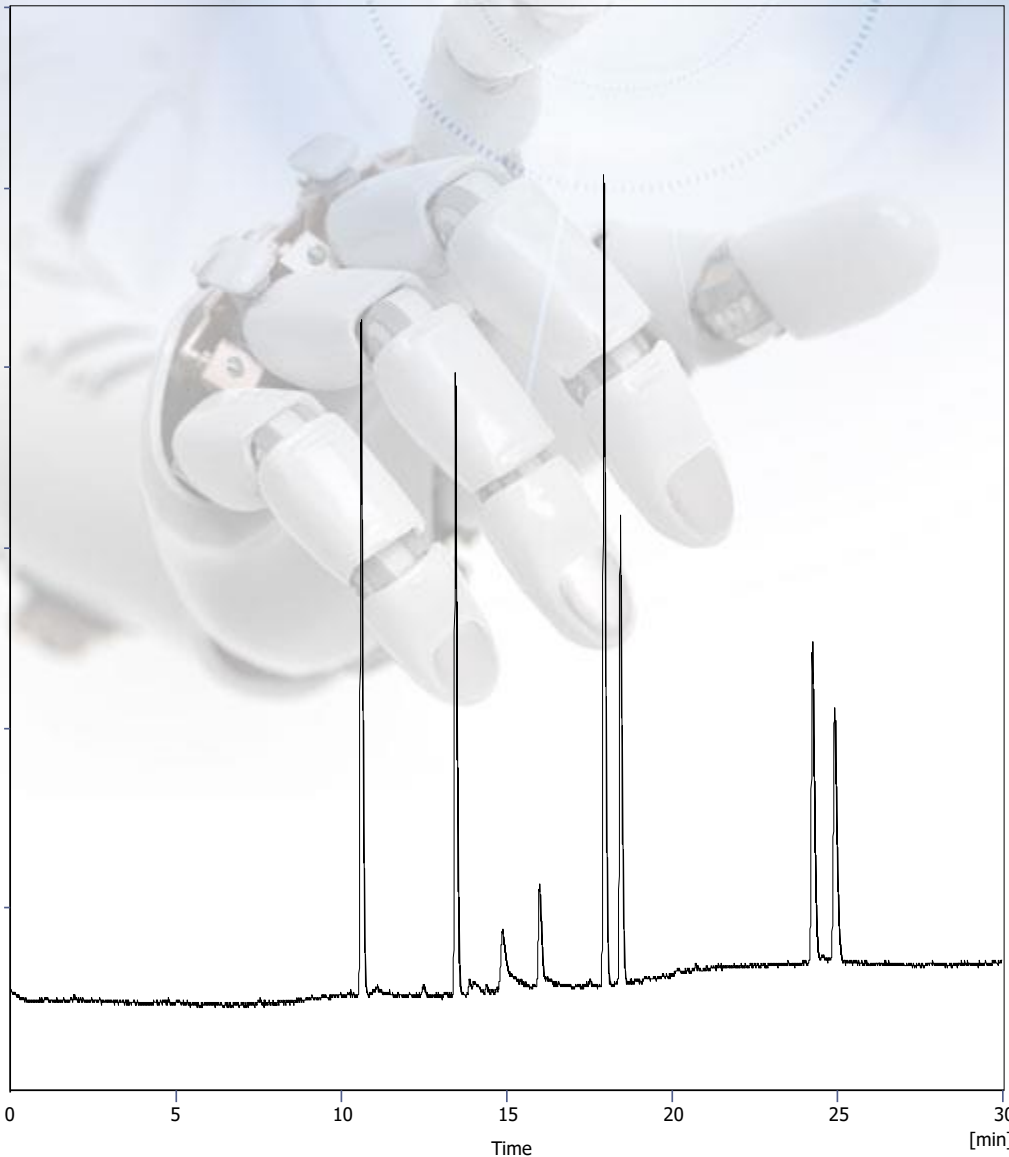