



Conditions

Column:	ACE Excel 5 C18
Dimensions:	250 x 4.6 mm
Part Number:	EXL-121-2546U
Mobile Phase:	MeOH/ H ₂ O/ 0.2 M NaH ₂ PO ₄ in H ₂ O/ 0.4 M tetrabutylammonium hydroxide in H ₂ O (450:447:100:3 v/v/v/v). Adjusted to pH 5.50 with H ₃ PO ₄
Flow Rate:	1.0 mL/min
Injection:	10 µL
Temperature:	40 °C
Detection:	UV, 254 nm
Monograph Method:	ChP.2015, Vol. 2

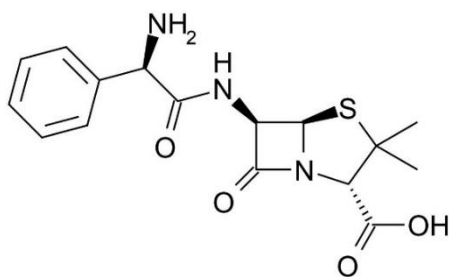

1. Ampicillin

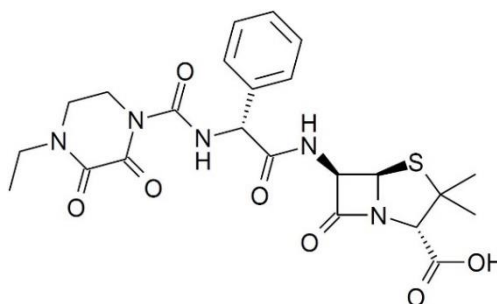

2. Piperacillin

