

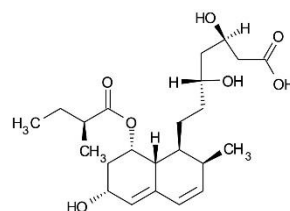


Conditions


Column: ACE 3 C18
Dimensions: 50 x 3.0 mm
Part Number: ACE-111-0503
Mobile Phase: MeCN/MeOH/THF/H₂O/Acetic Acid (15:20:5:60:0.1 v/v/v/v/v)
Flow Rate: 0.6 mL/min
Injection: 2 µL
Temperature: Ambient
Detection: API 3000 triple quad MS
TurbolonSpray – negative mode
Extracted ion chromatogram of MRM m/z 423.3 → 321.1

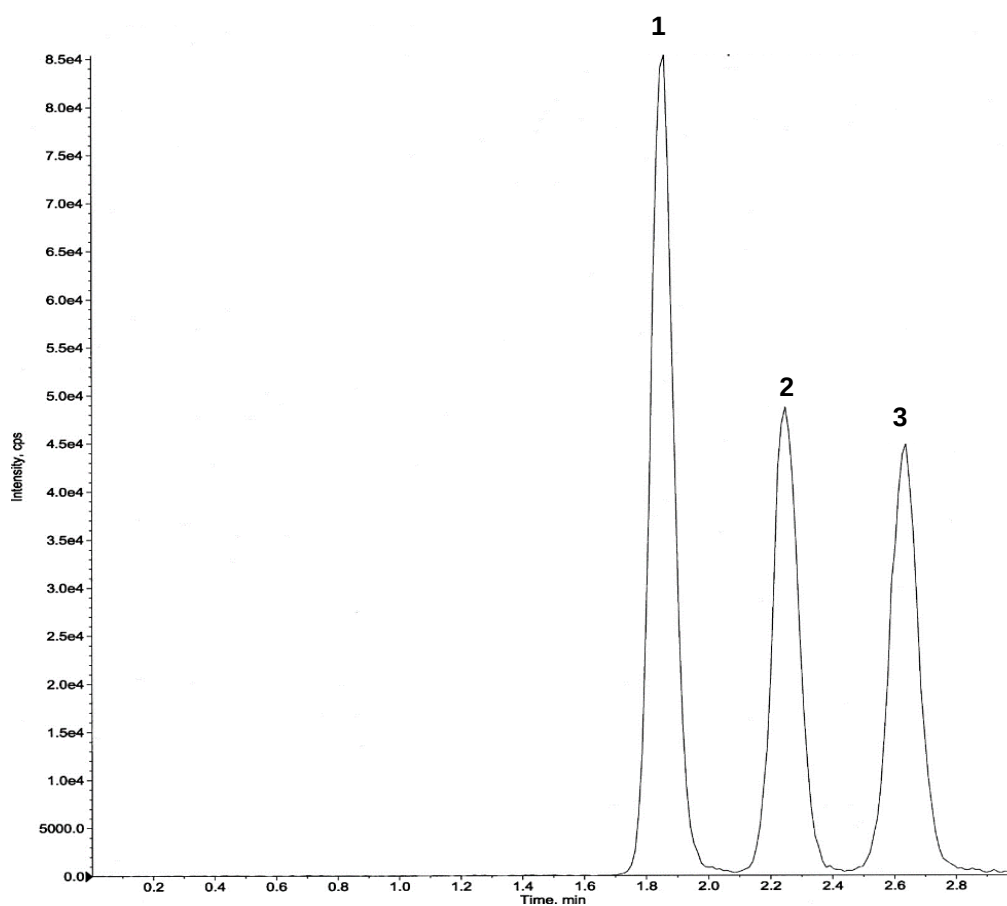

1. 6-epi-Pravastatin
(MW: 424)

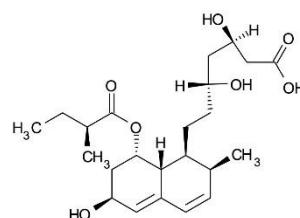

2. Pravastatin
(MW: 424)

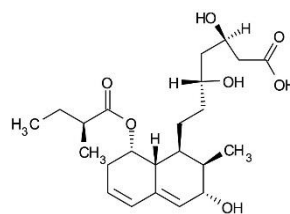

3. 3 α -iso-Pravastatin
(MW: 424)