
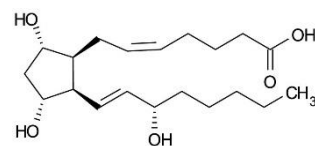


Conditions


Column: ACE 3 C18
Dimensions: 50 x 2.1 mm
Part Number: ACE-111-0502
Mobile Phase: A: H₂O
B: MeOH
C: MeCN

Time (mins)	%A	%B	%C
0	70.0	20.0	10.0
9	10.0	60.0	30.0
10	0.1	66.6	33.3

Flow Rate: 0.2 mL/min
Injection: 10 µL
Temperature: 40 °C
Detection: ESI (-) MS/MS
MRM 353.3 → 193

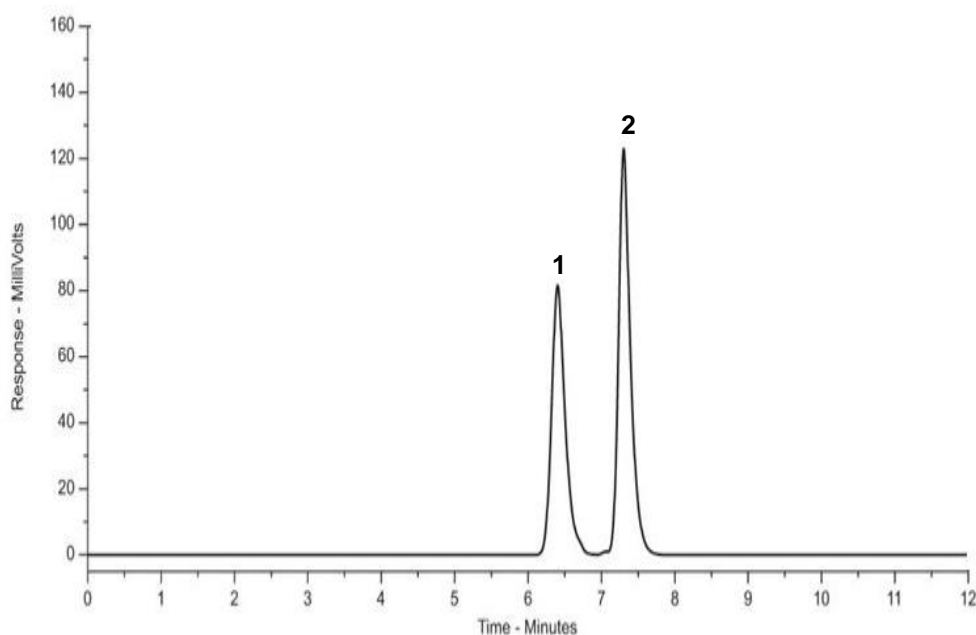

1. 8-iso-Prostaglandin F_{2α}

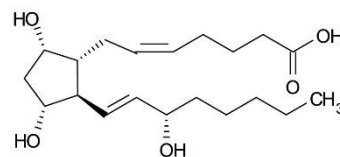

2. Prostaglandin F_{2α}